

qPCR guide



Eurogentec
Experience true partnership

PREFACE

Real-Time quantitative PCR, or qPCR in short, is heralded as the gold standard for accurate, sensitive and fast quantification of nucleic acid sequences. Indeed it is a wonderful technology, simple to perform and the risk for cross contamination is extremely low, as reaction tubes do not need to be opened. However, the simplicity of generating data is also the Achilles heel of qPCR. Nowadays, high quality reagents are available, and dedicated instruments do a perfect job. Hence, anyone that masters the skill of pipetting can generate beautiful sigmoidal amplification curves. What many investigators unfortunately fail to appreciate is quality control, which is essential throughout the entire qPCR workflow (from living cells, over extraction of nucleic acids, storage, various enzymatic steps such as DNase treatment, reverse transcription and PCR amplification, and finally data-analysis) in order to draw biologically meaningful conclusions. It is now well over 10 years after the conception of qPCR and there is still no consensus on good laboratory practice in general, and the level of quality control, experiment design and the data-analysis strategy in particular. As such, the life science literature is flawed with studies that are meaningless, inconclusive, or in the worst case erroneous (resulting in retraction of some papers). Fortunately, standards for data exchange and reporting guidelines are finding their way to the scientific community (see for example efforts made by the RDML consortium, <http://www.rdml.org>).



When inspecting the scientific literature relating to qPCR, it is clear that the number of publications over the years follows the same exponential rate as the PCR reaction itself. Currently, there is no plateau phase in sight; qPCR remains a growing market, especially in clinical diagnostics. In the research field, gene expression analysis is by far the most popular application, with double stranded DNA specific binding dyes and hydrolysis probes as the dominating detection chemistries (source: the Real-Time PCR Primer and Probe database, RTPPrimerDB, <http://medgen.ugent.be/rtprimerdb/>).

Coming from relatively low throughput qPCR systems a decade ago, 384-well instruments have now become more and more mainstream. “New kids” are now also on the block, implementing nanofluidic technologies, significantly increasing throughput and decreasing reaction volumes (down to nanoliter range), without the need for advanced liquid handling instruments. Following this promising path of ever increasing throughput, we might foresee a future in which we will simply PCR profile entire signalling networks or pathways, or even determine the expression level of almost every gene in an organism within only a few hours, with unprecedented sensitivity, accuracy and at low cost. Surely, it’s a long road with many hidden obstacles, but witnessing the tremendous progress the science community have made in the last decade, the envisioned whole genome PCR profiling comes within reach.

Jo Vandesompele
 Professor functional genomics and applied bio-informatics, Ghent University
 Co-founder and CEO, Biogazelle, the Real-Time PCR data analysis company



TABLE OF CONTENTS

Introduction	6
Preliminary steps	7
Template preparation: DNA or RNA	7
One-step or two-step reaction?	8
Which type of primers for the Reverse Transcription?	9
Qualitative vs quantitative PCR technologies	10
What are qualitative PCR & quantitative PCR?	10
Quantitative PCR	11
Qualitative PCR	11
Detection methods	12
SYBR® Green I	12
High Resolution Melting dyes (HRM dyes)	13
TaqMan® probes = Double-Dyes probes	13
LNA® Double-Dye probes	15
Molecular Beacon probes	16
Scorpions® primers	18
Hybridization probes (also called FRET probes)	19
TaqMan® MGB® probes	20
MGB Eclipse® probes	20
Other technologies	21
When should I multiplex?	21
How do you get started in qPCR?	22
Thermocycler	22
qPCR reagents : formats and components?	24
Core kits or MasterMixes?	24
Role of core kits and MasterMixes components	24
Probe and primer quality	26
Manual coupling	26
Automatic coupling	26
Purification level	26

Guidelines for primers and probes design	27
General rules	27
Design guidelines for SYBR® Green I assays	28
Design guidelines for Double-Dye Probe assays	29
Quenching principle	31
Fluorophore and quencher: definitions	31
Different ways of quenching	31
FRET quenching	32
Optimal fluorophore quencher combination	32
Dye-thermocyclers compatibility table	34
Double-Dye Oligonucleotides in single-color detection	34
Double-Dye Oligonucleotides for multiplexing detection	34
Example of a classical assay	36
Some tips and tricks before starting	36
Classical protocol and thermal profile	37
Useful definitions	38
Normalization and quantification methods	39
Absolute quantification	41
Relative quantification	41
Evaluation of your qPCR assay	43
Controls	43
Negative controls	43
Positive controls	44
Commonly used control kits	45
Singleplex or multiplex assay for the controls?	45
Role of the standard curve and evaluation of its parameters	46
PCR efficiency	46
R square (R^2)	47
Sensitivity	47
Reproducibility	47
Specificity of a SYBR® Green I assay	47
Optimization tips	49
In general	49
For probe multiplexing assay/qualitative qPCR	50
Troubleshooting guide	51
Frequently Asked Questions	53
Eurogentec products	55
References	64



This guide is dedicated for both beginners and experts in qPCR. It explains both the basics of qPCR and gives useful tips for troubleshooting. Real-Time qPCR is a technique that involves a high number of parameters, all of them having to be optimized to get the highest quality assay, and all of which are specific to one assay only. There is consequently no “magic formula”. However, this booklet can provide tools to help you, from your first contact with qPCR to the optimization steps. If you are an expert, just go straight to the last pages of this booklet. We hope that by working through the pages, we will help you enjoy this technique and reach the best results for your research!

INTRODUCTION

The Polymerase Chain Reaction (PCR) has been invented in 1983 by Kary Mullis (Nobel Prize in 1993), (Mullis K. et Fobona F., 1987). Three years after its invention, there was an incredible expansion of its use thanks to the commercialization of the *Taq* polymerase, a polymerase that resists high temperatures. In 1991, the first hydrolysis probe was used in combination with the technique. In 1992, the technique was again improved by the use of Ethidium Bromide (EtBr), thanks to the fluorescence that results from the binding to duplex DNA. The kinetics of fluorescence accumulation during thermocycling was directly related to the starting number of DNA copies. This was the starting point of Real-Time qPCR. Today, 25 years after Mullis’s discovery, both PCR and qPCR are widely used technologies. The principle, and aim, of the PCR technology is to specifically increase a target from an undetectable amount of starting material. In classical PCR, at the end of the amplification, the product can be run on a gel for detection of this specific product. In Real-Time PCR, this step can be avoided since the technology combines the DNA amplification with the immediate detection of the products in a single tube. The homogeneous format is highly beneficial as it removes the significant contamination risk caused by opening tubes for post-PCR manipulation. It is also less time consuming than gel based analysis and can give a quantitative result (Kubista M. *et al.*, 2006).

Current detection methods are based on changes in fluorescence, which are proportional to the increase of target. Fluorescence is monitored during each PCR cycle providing an amplification plot, allowing the user to follow the reaction in real time. (Figure 1).

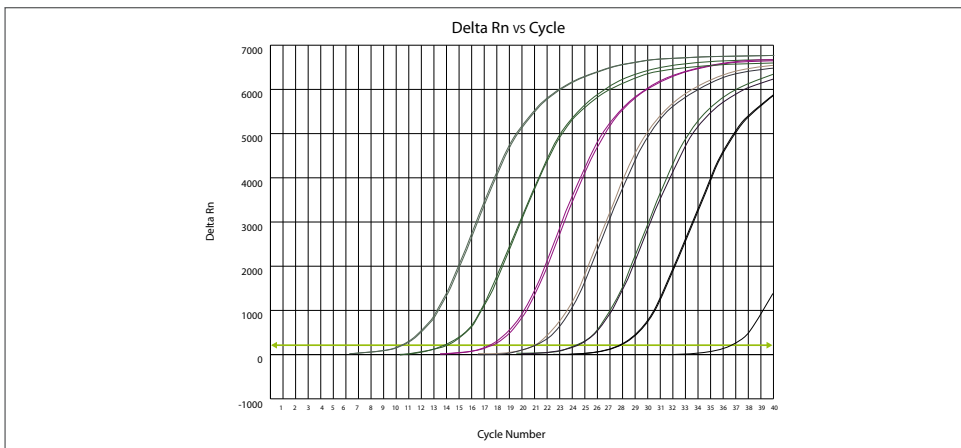


Figure 1. Amplification plot of 10 x serial diluted cDNA, using SYBR® Green detection method, on an ABI 7500 thermocycler.

PRELIMINARY STEPS

Variability in qPCR is often related to steps upstream to the qPCR step itself. Sample extraction, quantity of sample, or efficiency of the reverse transcription (RT) are some of the many parameters that can influence the results of your qPCR assay.

That's why it is important to also consider the quality of these steps prior to performing your qPCR assay (Fleige S. and Pfaffl MW, 2006).

Template preparation: DNA or RNA

DNA and RNA can be extracted with Trizol. This includes many steps and a careful wash of the sample, hence we recommend the use of an appropriate commercially available kit. Many column-based kits, which will contain all the required reagents for the full extraction / purification procedure are available.

These kits will also outline general guidelines, such as the storage conditions and shelf life of the extracted RNA or DNA. However these guidelines may vary between kits due to the different composition of buffers. It is also recommended that the buffers supplied with each kit are used according to the manufacturer protocol.

The final product should be cleaned and free from any residual buffers such as EDTA or buffers containing solvents. Otherwise, these residual products may inhibit the action of the *Taq* polymerase for the PCR step, or could modify the salt concentration of the buffer. This is usually not a problem when using spin column kits instead of manual extraction techniques.

It is worth bearing in mind, that this first step is nearly the most important of your RT-qPCR assay, as the quality of the extraction will influence the quality of your detection and quantification. It's mostly important to ensure the reproducibility of your extraction steps, if you wish to compare biological samples.

The main problems that can occur, are extraction of inhibitors together with the nucleic acids or degradation of the sample. Inhibitors are mainly found in blood sample or environmental samples. For example, humic acid can be a strong inhibitor in samples extracted from soil. As a general rule, it is always better to extract DNA/RNA from fresh sample and to store it at -80 °C.



One-step or two-step reaction?

If your starting material is total RNA, you have to perform the Reverse Transcription (RT) step before your qPCR assay. This can be done in One-Step or Two-Steps. Both methods have pros and cons you might consider to choose the best one for your application.

In both cases, the RT reaction should be set up in a clean environment to avoid contamination. We also recommend working with RNase free plastics. The tubes containing the reaction should be maintained on ice during the set up of the reaction. This will prevent an early reaction, generated before the incubation period.

The RT step lasts 30 minutes around 50 °C. However, if random primers or oligo d(T) are used (two-step reaction), an initial 10 min-step at 25 °C is necessary to maximize the annealing efficiency, as the short oligo d(T) have a lower T_m. The next paragraph explains what kind of primers can be used under which circumstances.

The cDNA generated should be stored at -80 °C. The use of a separate RT step is recommended when the reaction is performed with a limiting amount of starting material.

	One-Step qRT-PCR	Two-step qRT-PCR
Description	<p>RT and qPCR reactions performed in the same tube</p> <p>Optimized working buffer for both the RT and PCR enzymes</p>	<p>RT and qPCR reactions performed in separate tubes</p>
Pros	<p>Save pipetting steps</p> <p>No contamination between RT and qPCR steps</p> <p>Lower background in SYBR® assays</p> <p>Best option for High-Throughput Screening (less time consuming than 2-step reactions)</p>	<p>More efficient because random primers and oligo d(T) can be used</p> <p>Possibility to stock cDNA to quantify several targets</p> <p>More flexible (separate optimization possible for the two reactions)</p>
Cons	<p>No possibility to use UNG carry-over prevention</p> <p>Usually less sensitive than a Two-Steps qRT-PCR assay</p>	<p>RNase inhibitors that can influence the PCR reaction after the RT</p> <p>Higher background when performing a SYBR® assay</p>

Table 1: Pros and cons of One-Step qRT-PCR versus Two-Step qRT-PCR

Which type of primers for the Reverse Transcription?

When performing a RT, it is possible to use three different types of primers:

- short random primers, which bind anywhere in the genome and allow the reverse transcriptase to fill up the gaps, lead to higher yields
- oligo d(T) bind to the poly-A tail of the RNA and then only transcribe RNA. This will avoid contamination by genomic DNA. As the poly-A tail is located at the extremity of the gene, it will also lead to more full transcripts
- specific primers, which bind to the gene of interest

The combination of oligo d(T) primers and random primers will give the highest yields and the longest transcripts, whereas specific primers transcribe only specific RNA but reduce the yield.

With a One-Step qRT-PCR reaction, using specific primers is the only option, as it should be avoided that oligo d(T) primers and random nonamers participate in the PCR step, giving many aspecific products. The RT step is performed at 40-50 °C and the specific PCR primers are designed for the PCR step at 60 °C, leading to partial annealing of the primers during the RT step. In a Two-Steps qRT-PCR kit, oligo d(T) primers or/and random nonamers are used in the RT step, and specific primers in the PCR step, leading to specific cDNA.



QUALITATIVE vs. QUANTITATIVE PCR

What are qualitative PCR & quantitative PCR?

Firstly it is important to review the principles of PCR, before making the distinction between qualitative and quantitative PCR.

The action of the *Taq* polymerase allows extension of short single-stranded synthetic oligonucleotides (i.e. primers) during repeated cycles of heat de-naturation, primer annealing and primer extension.

The primers are designed to specifically bind the extremities of the DNA fragment to be amplified. The *Taq* polymerase uses the target DNA added to the reaction as a template for primer extension. At each cycle, more DNA is synthesized, creating more template DNA. The reaction proceeds in an exponential manner, doubling the amount of target during each cycle, until one of the reagents becomes limiting and the reaction reaches a plateau.

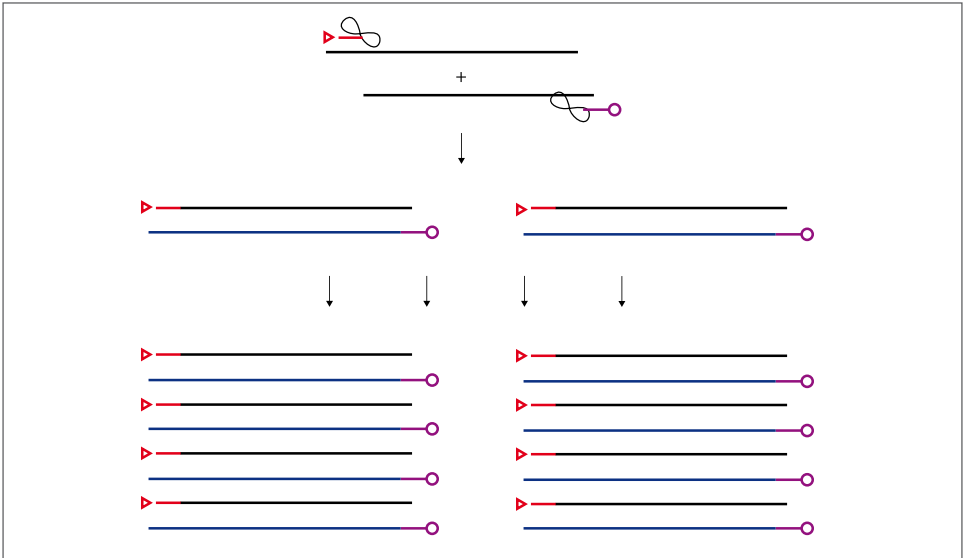


Figure 2. The Polymerase Chain Reaction. (K. Mullis and F. Fabona, 1987). The red and pink lines represent the primers that the polymerase extends from. The blue and black lines represent the single-stranded DNA template produced from de-naturation of the double-stranded DNA template.

Quantitative PCR

In quantitative qPCR, a specific or non-specific detection chemistry allows the quantification of the amplified product. The amount detected at a certain point of the run is directly related to the initial amount of target in the sample.

The most common applications of quantitative PCR are gene expression analysis, pathogen detection/quantification and microRNA quantification (Schmittgen TD *et al.*,2008).

For example, effects of different treatments, on the level of mRNA transcription can be measured. Quantitative PCR software uses the exponential phase of PCR for quantification. PCR is initially an exponential process but eventually reaches a plateau phase, when one of the reagents becomes limited. Reactions can plateau at different levels even if they have the same starting concentration of target. During the exponential phase, the amount of target is assumed to be doubling every cycle and no bias is expected due to limiting reagents. Analysis takes the Ct (cycle number) value, at the point when the signal is detected above the background and the amplification is in exponential phase. The more abundant the template sample, the quicker this point is reached, thus giving earlier Ct values. Differences in Ct then have to be correlated to some other quantitative values to make them meaningful (*cf.* p. 39).

Qualitative PCR

In qualitative qPCR, the goal is to detect the presence or absence of a certain sequence. It could be for virus sub-typing and bacterial species identification for example. It can also be used for allelic discrimination between wild type and mutant, between different SNPs (Single Nucleotide Polymorphisms) or between different splicing forms. In this case, different fluorophores can be used for the two alleles, and the ratio of the fluorophores signals correlates to the related amount of one form compared to the other one.

Specific detection methods such as Double-Dye probe systems are more often used for these applications, and probes can be used to detect single base mutations or small deletions.

Recently, a new technique called HRM (High Resolution Melting) emerged, allowing gene scanning and SNP genotyping without using probes (*cf.* p. 13).

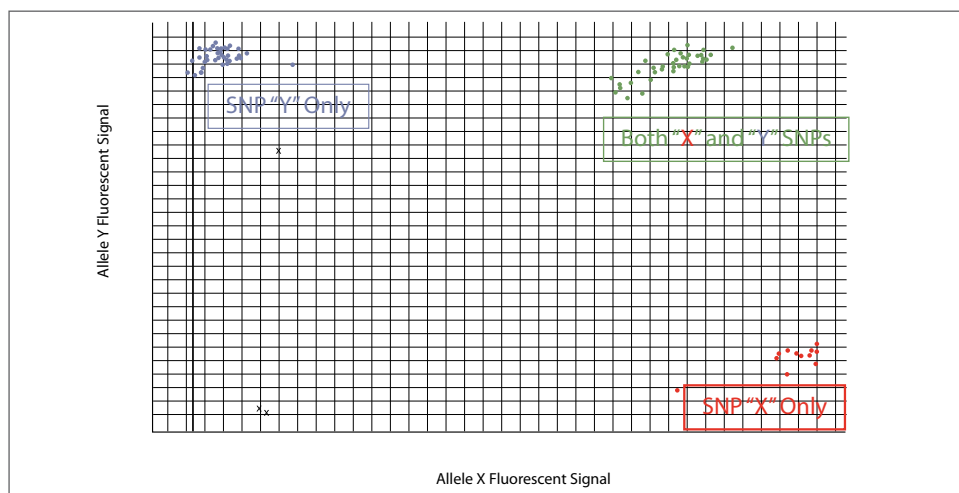


Figure 3. Allelic discrimination using two Double-Dyes labelled with different fluorophore (*ncvs.org*).



Detection methods

In order to detect and measure the amount of target in the sample, a measurable signal has to be generated, which is proportional to the amount of amplified product. All current detection systems use fluorescent technologies.

Some of them are non-specific techniques, and consequently only allow the detection of one target at a time. Alternatively, specific detection chemistries can distinguish between non-specific amplification and target amplification. These specific techniques can be used to multiplex the assay, i.e. detecting several different targets in the same assay.

SYBR® Green I

SYBR® Green I is the most commonly used dye for non-specific detection. It is a double-stranded DNA intercalating dye, that fluoresces once bound to the DNA. A pair of specific primers is required to amplify the target with this chemistry. The amount of dye incorporated is proportional to the amount of generated target. The dye emits at 520 nm and fluorescence emitted can be detected and related to the amount of target.

The inconvenience of this technique is that the SYBR® Green I will bind to any amplified dsDNA. Consequently, primer dimers or unspecific products introduce a bias in the quantification. However, it is still possible to check for the specificity of the system by running a meltcurve at the end of the PCR run (*cf.* p. 47). The principle is that every product has a different dissociation temperature, depending of the size and base contents, so it is still possible to check the number of products amplified. A valid SYBR® assay - primer pair - should produce a unique, well defined peak on the meltcurve.

For these reasons, SYBR® Green I is rarely used for qualitative PCR. However, SYBR® Green I is often used as the first step to optimize a specific detection system assay, to check the specificity of the primers and validate the design.

SYBR® Green I advantages

- Low cost assay
- Easy design and set up

SYBR® Green I disadvantages

- Non specific system
- Not adapted to multiplex
- Non suitable for qualitative qPCR

High Resolution Melting dyes (HRM dyes)

High Resolution Meltcurve analysis is a newly emerging technology, which characterizes nucleic acid samples based on their dissociation behaviour. It combines the principle of intercalating dyes, meltcurve analyses and the application of specific statistical analyses.

HRM uses the fundamental property of the separation of the two strands of DNA with heat (melting), and the monitoring of this melting with a fluorescent dye. On the contrary of SYBR Green, HRM dyes do not inhibit PCR at high concentration. The dye can consequently saturate the amplified target dsDNA and fluoresces. Melting temperature of a dsDNA target depends on GC content, length, and sequence. Due to the high sensitivity of HRM dyes, even a single base change will induce differences in the melting curve, and consequently in fluorescence (Erali M. *et al.*, 2008).

Main applications of HRM include gene scanning (search for the presence of unknown variations in PCR amplicons), SNP genotyping, DNA methylation analysis, DNA mapping, DNA fingerprinting.

This emerging method is less expensive and as precise than probe-based methods.

Only a few thermocyclers on the market currently allow the use of this technology, among them the Roche LightCycler®480, the Corbett Life Science Rotor-Gene™ 6000, and the ABI Prism® 7500. The main HRM dyes available are EvaGreen, LCGreen®, SYTO® 9 and BEBO.

TaqMan® probes = Double-Dye probes

TaqMan® probes, also called Double-Dye Oligonucleotides, Double-Dye Probes, or Dual-Labelled probes, are the most widely used type of probes and are often the method of choice for scientists who have just started using Real-Time PCR. They were developed by Roche (Basel, Switzerland) and ABI (Foster City, USA) from an assay that originally used a radio-labelled probe (Holland *et al.* 1991), which consisted of a single-stranded probe sequence that was complementary to one of the strands of the amplicon.

A fluorophore is attached to the 5' end of the probe and a quencher to the 3' end. The fluorophore is excited by the machine and passes its energy, via FRET (Fluorescence Resonance Energy Transfer) to the quencher. Traditionally the FRET pair has been FAM as the fluorophore and TAMRA as the quencher. In a well designed probe (*cf.* Taqman® design rules on p. 30), FAM does not fluoresce as it passes its energy onto TAMRA. As TAMRA fluorescence is detected at a different wavelength to FAM, the background level of FAM is low. The probe binds to the amplicon during each annealing step of the PCR. When the *Taq* polymerase extends from the primer which is bound to the amplicon, it displaces the 5' end of the probe, which is then degraded by the 5'-3' exonuclease activity of the *Taq* polymerase. Cleavage continues until the remaining probe melts off the amplicon (*cf.* fig. 4).

This process releases the fluorophore and quencher into solution, spatially separating them (compared to when they were held together by the probe). This leads to an irreversible increase in fluorescence from the FAM and a decrease in the TAMRA.



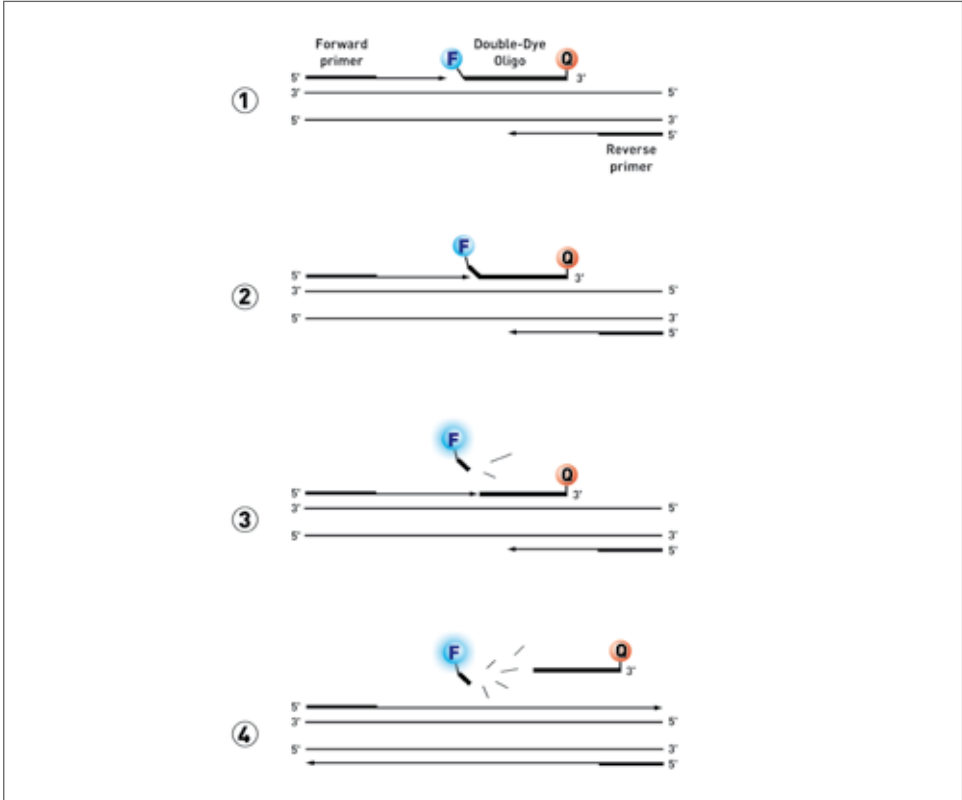


Figure 4. TaqMan® probes mode of action

TaqMan® probes can be used for both quantification and mutation detection, and most designs appear to work well. They are convenient either for allelic discrimination or expression profiling and are usually easy to design, easy to use in a standard protocol and need minimal optimization. For mutation detection, the probe is designed to hybridize over the mutation site and can be made specific enough to detect single base differences.

To obtain robust allelic data, it is vital that a different probe is used for each different mutation, otherwise negative results with one probe, (caused by failed PCR reactions), can be scored as an absence of a particular allele. Ideally these probes would be multiplexed to make the assays faster and cheaper.

Double-Dye probes advantages

- Widely used, several modifications possible
- Specific system
- Multiplex capabilities

Double-Dye probes disadvantages

- More expensive than SYBR® Green I, but the lowest cost specific system available

LNA® Double-Dye probes

LNA® (Locked Nucleic Acid) was developed by Exiqon® (Vedbaek, Denmark). LNA® changes the conformation of the helix and increases the stability of the duplex. The integration of LNA® bases into Double-Dye Oligonucleotide probes, opens up great opportunities to improve techniques requiring high affinity probes as specific as possible, like SNP detection, expression profiling and *in situ* hybridization.

LNA® is a bicyclic RNA analogue, in which the ribose moiety in the sugar-phosphate backbone is structurally constrained by a methylene bridge between the 2'-oxygen and the 4'-carbon atoms.

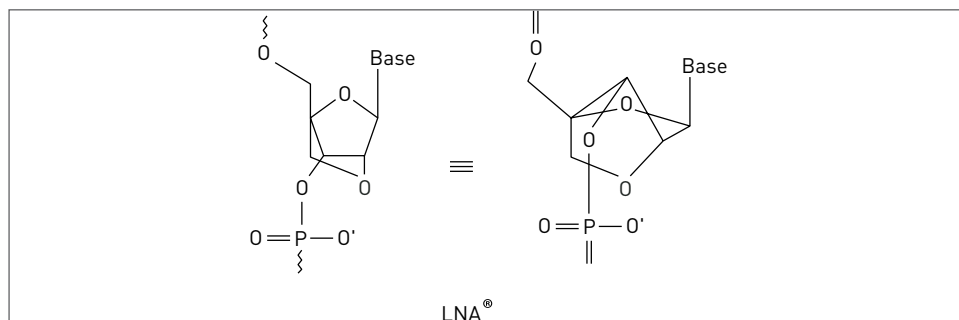


Figure 5. LNA® backbone structure

The integration of LNA® bases into probes changes the conformation of the double helix from the B to A type (Ivanova A. *et al.*, 2007).

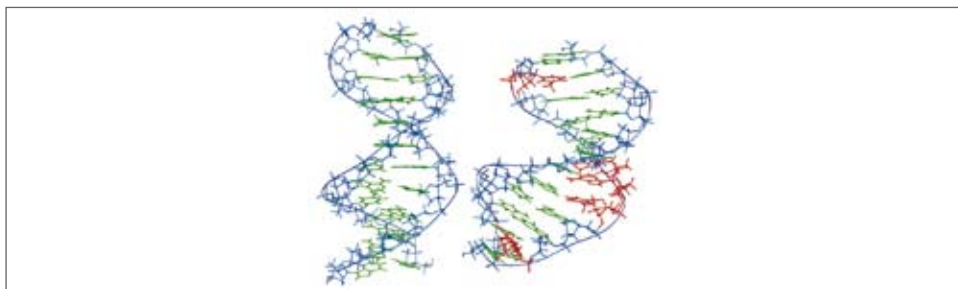


Figure 6. B helix and A helix containing LNA® bases

LNA® conformation allows a much better stacking and therefore a higher stability. By increasing the stability of the duplex, the integration of LNA® monomers into the oligonucleotide sequence allows an increase of the melting Temperature (T_m) of the duplex. It is therefore possible to reduce the size of the probe, which increases the specificity of the probe and helps designing it (Karkare S. *et al.*, 2006).

LNA® advantages

- Increases thermal stability towards complementary DNA and RNA
- Allows use of shorter probes
- High specificity and reproducibility
- Simplifies multiplex assays by adjusting T_m values of primers and probes



Molecular Beacon probes

Molecular Beacons are probes that contain a stem-loop structure, with a fluorophore and a quencher at their 5' and 3' ends, respectively. The stem is usually 6 bases long, should mainly consist of C's and G's, and holds the probe in the hairpin configuration (Li Y. *et al.*, 2008).

The 'stem' sequence keeps the fluorophore and the quencher in close vicinity, but only in the absence of a sequence complementary to the 'loop' sequence. As long as the fluorophore and the quencher are in close proximity, the quencher absorbs any photons emitted by the fluorophore. This phenomenon is called collisional (or proximal) quenching. In the presence of a complementary sequence, the Beacon unfolds and hybridizes to the target, the fluorophore is then displaced from the quencher, so that it can no longer absorb the photons emitted by the fluorophore, and the probe starts to fluoresce. The amount of signal is proportional to the amount of target sequence, and is measured in real time to allow quantification of the amount of target sequence (Takacs T. *et al.*, 2008).

The increase in fluorescence that occurs is reversible, (unlike TaqMan® probes), as there is no cleavage of the probe, that can close back into the hairpin structure at low temperature.

The stem structure adds specificity to this type of probe, because the hybrid formed between the probe and target has to be stronger than the intramolecular stem association. Good design of Molecular Beacons can give good results, however the signal can be poor, as no physical separation of fluorophore from quencher occurs.

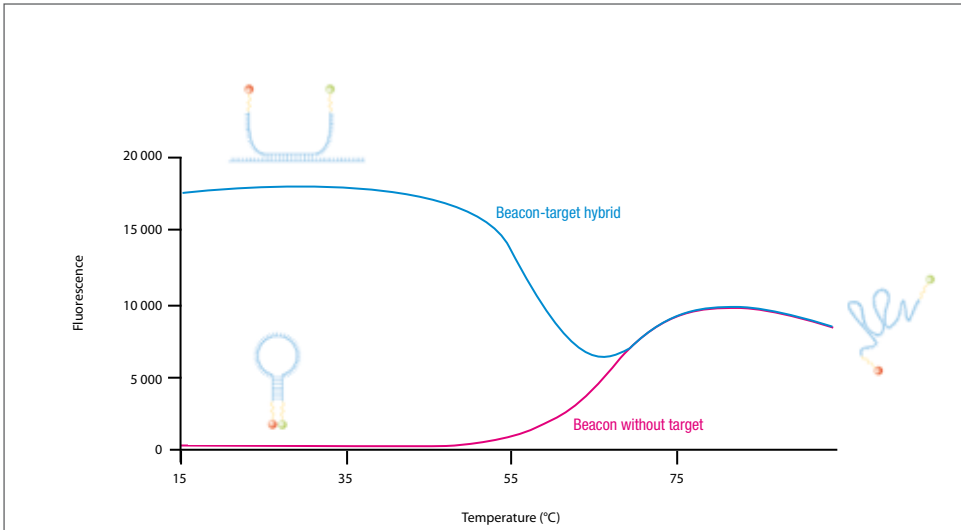


Figure 7. Melt curve of a Molecular Beacon with and without a synthetic complement as a target. Annotations show the configuration of the beacon and target at different temperatures.

Two improved forms of Molecular Beacons exist:

- Wavelength-Shifting Molecular Beacons**

Wavelength-Shifting Molecular Beacons are brighter than standard Molecular Beacons due to an enhanced fluorescence intensity of the emitter fluorophore. These probes contain a harvester fluorophore that absorbs strongly in the wavelength range of the monochromatic light source, an emitter fluorophore of the desired emission color, and a non-fluorescent (dark) quencher. In the absence of complementary nucleic acid targets, the probes are non-fluorescent, whereas in the presence of targets, they fluoresce, not in the emission range of the harvester fluorophore, that absorbs the light, but rather in the emission range of the emitter fluorophore. This shift in emission spectrum is due to the transfer of the absorbed energy from the harvester fluorophore to the emitter fluorophore by FRET, which only takes place in probes that are bound to the targets. Wavelength-Shifting Molecular Beacons are substantially brighter than conventional Molecular Beacons that cannot efficiently absorb energy from the available monochromatic light source (Tyagi S. *et al.*, 2000).
- 2' O-Methyl RNA Molecular Beacons**

To detect the various RNA classes in living cells, several approaches have been developed. One of these is based on the use of 2'-O-Methyl RNA Molecular Beacon probes for the detection of small nuclear RNAs. 2'-O-methyl RNA probes are considered to perform better than DNA oligonucleotides because they are not only nuclease resistant, but also possess a higher affinity, increased specificity, faster hybridization kinetics, and a superior ability to bind to structured targets compared to DNA oligonucleotides. Recently Molecular Beacons were introduced for RNA detection in living cells. The rationale for using Molecular Beacons to detect RNAs in living cells was to improve signal-to-noise ratios by eliminating fluorescence signals derived from non-hybridized probe sequences. It appears that linear 2'-O-Methyl RNA probes are usually more suitable for specific detection of these RNAs, representing different classes of RNA, in the nuclei of living cells. Molecular Beacons result in images with improved signal to noise ratios, thereby leading to better detection sensitivity.

Molecular Beacons

- Reversible phenomenon (melting curve possible)
- Increased discriminatory competence due to competition intra-inter target



Scorpions[®] primers

Scorpions[®] primers are suitable for both quantitative Real-Time PCR and genotyping/end-point analysis of specific DNA targets.

They are PCR primers with a “stem-loop” tail consisting of a specific probe sequence, a fluorophore and a quencher. The “stem-loop” tail is separated from the PCR primer sequence by a “PCR blocker”, a chemical modification that prevents the *Taq* polymerase from copying the stem loop sequence of the Scorpions[®] primer. Such read-through would lead to non-specific opening of the loop, causing a non-specific fluorescent signal.

The hairpin loop is linked to the 5' end of a primer via a PCR blocker. After extension of the primer during PCR amplification, the specific probe sequence is able to bind to its complement within the same strand of DNA. This hybridization event opens the hairpin loop so that fluorescence is no longer quenched and an increase in signal is observed.

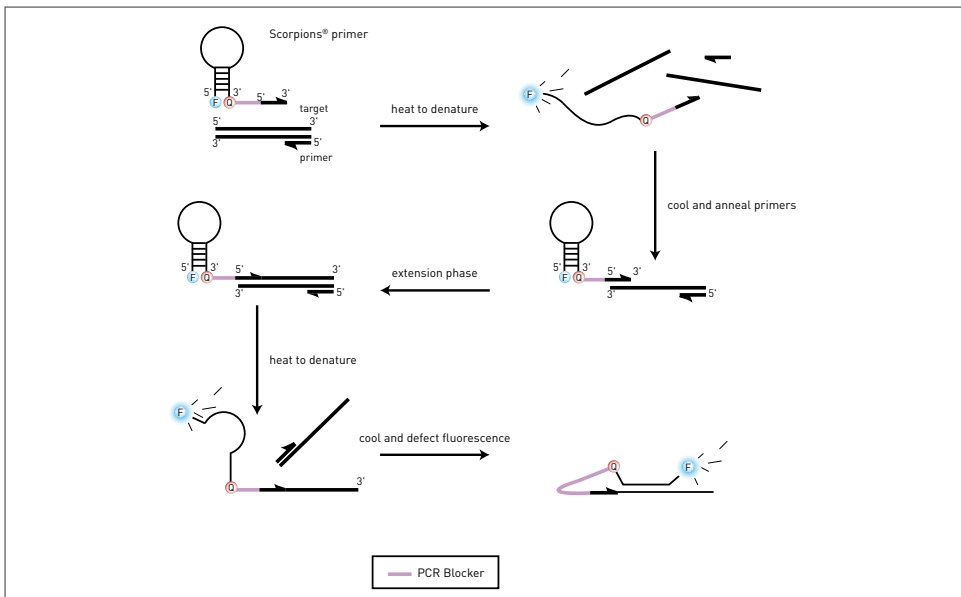


Figure 8. Scorpions[®] unimolecular mode of action. Progression through a PCR reaction.

Unimolecular probing is kinetically favorable and highly efficient. Covalent attachment of the probe to the target amplicon ensures that each probe has a target in the near vicinity. Enzymatic cleavage is not required, thereby reducing the time needed for signaling compared to TaqMan[®] probes, which must bind and be cleaved before an increase in fluorescence is observed.

Scorpions[®] primers have successfully been used for mutation detection and quantification (Solinas A. *et al.*, 2001), having the specificity and melt curve analysis of the Molecular Beacons, with additional speed and efficiency. However, constraints associated with Scorpions[®] design is probably the reason why this type of probe is not as successful as it should be.

There are three types of Scorpions® primers :

- Standard Scorpions®, which consist of a bi-labelled probe with a fluorescent dye at the 5' end and an internal non-fluorescent quencher.
- FRET Scorpions®, for use on a LightCycler® system. As the capillary system will only excite at 470 nm (FAM absorption wavelength) it is necessary to incorporate a FAM within the stem. A ROX is placed at the 5' end of the Scorpions® primer, FAM is excited and passes its energy onto the ROX.
- Duplex Scorpions® have also been developed to give much better signal intensity than the normal Scorpions® format. In Standard Scorpions® the quencher and fluorophore remain within the same strand of DNA and some quenching can occur even in the open form. In the Duplex Scorpions® the quencher is on a different oligonucleotide and physical separation between the quencher and fluorophore is greatly increased, reducing the quenching when the probe is bound to the target.

Advantages of Scorpions® primers vs Taqman® probes

- Stronger signal and lower background
- Improved SNP applications
- Compatible with any dye, convenient for extended multiplexing assays
- Primer and probe on one molecule, kinetically favorable, fast cycling conditions

Hybridization probes (also called FRET probes)

Roche has developed hybridization probes (Caplin *et al.* 1999) for use with their LightCycler®. Two probes are designed to bind adjacent to one another on the amplicon. One has a 3' label of FAM, whilst the other has a 5' LC dye, LC red 640 or 705.

When the probes are not bound to the target sequence, the fluorescent signal from the reporter dye is not detected. However, when the probes hybridize to the target sequence during the PCR annealing step, the close proximity of the two fluorophores allows energy transfer from the donor to the acceptor dye, resulting in a fluorescent signal that is detected.

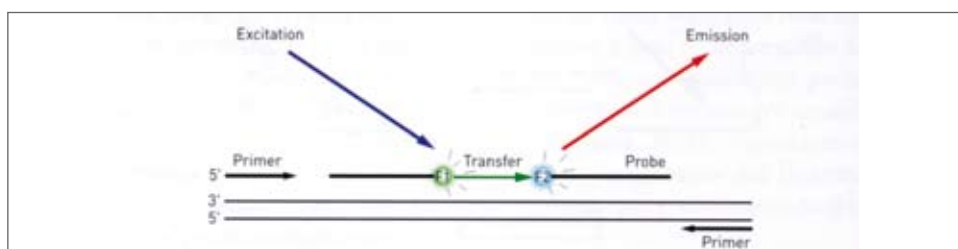


Figure 9. FRET probe principle



TaqMan® MGB® probes

TaqMan® MGB® probes have been developed by Epoch Biosciences (Bothell, USA) and Applied Biosystems (Foster City, USA). They bind to the minor groove of the DNA helix with strong specificity and affinity. When the TaqMan® MGB® probe is complemented with DNA, it forms a very stable duplex with DNA. The probe carries the MGB® moiety at the 3' end. The MGB strongly increases the probe T_m , allowing shorter, hence more specific designs.

The probe performs particularly well with A / T rich regions, and is very successful for SNP detection (Walburger *et al.*, 2001). It can also be a good alternative when trying to design a probe which should be located in the splice junction (for which conventional probes are hard to design). Smaller probes can be designed with T_m as 65-67 °C, which gives a better discrimination (the probe is more specific for single mismatch). A good alternative to MGB probes are LNA® probes where the increase in T_m induced by the addition of LNA® bases is specific, contrary to the MGB moiety (*cf.* p. 15).

During the primer extension step, the hybridized probe is cleaved by the 5' exonuclease activity of *Taq* polymerase and an increase in fluorescence is seen. Fluorescence of the cleaved probe during PCR is monitored in Real-Time by the thermocycler.

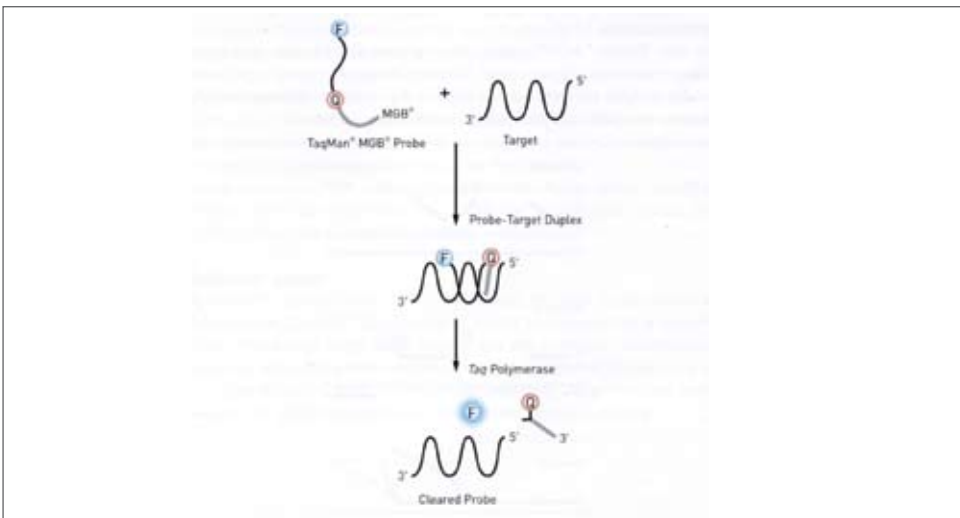


Figure 10. TaqMan® MGB® probe principle

MGB Eclipse® probes

MGB Eclipse® probes also known as QuantiProbes, have originally been developed by Epoch Biosciences (Bothell, USA). MGB Eclipse® probes carry a minor groove binder moiety that allows the use of short probes for very high specificity. These are short linear probes that have a minor groove binder and a quencher on the 5' end and a fluorophore on the 3' end. This is the opposite orientation to TaqMan® MGB® probes and it is thought that the minor groove binder prevents the exonuclease activity of the *Taq* polymerase from cleaving the probe. The quencher is a Non Fluorescent Quencher also known as Eclipse Dark Quencher. Quenching occurs when the random coiling of the probe in the free form brings the quencher and the fluorophore close to another. The probe is straightened out when bound to its target and quenching is decreased, leading to an increase in fluorescent signal.

Other technologies

The technologies that have been discussed above are the most widely used today, but numerous other technologies have occurred in publications, or are available on the market, such as: Resonsense probes, Light-up probes, HyBeacon[®] probes, LUX primers, Yin-yang probes, or Amplifluor[®]. You can contact us for more information on any of them.

When should I multiplex?

We have previously seen that some chemistries are more adapted to multiplexing than others. Theoretically, multiplex assay can be performed also with a non-specific system such as SYBR[®] Green I assay. The most common use of multiplex reaction is for qualitative PCR assay. In a quantitative assay, it mainly allows to include a control gene in the same well as the target (*cf.* p. 39, Normalization and Quantification methods).

Multiplexing is sometimes a necessity (SNP detection), but multiplex assays can also be set up for practical reasons, to save time and reagents (Ishii T. *et al.*, 2007).

When is it worth multiplexing? When is it better to singleplex?

- For SNP/mutation analysis: for bi-allelic discrimination, one specific probe is designed for each allele. Both of them have a 5' end fluorophore and a quencher at the 3' end and duplex assay is necessary.
- When targeting two genes with equal expression levels with PCR reactions of similar amplification efficiency, multiplex Real-Time PCR can be performed without any doubt.
- When analyzing two genes with different expression levels we recommend separating amplification of targets.
- If you quantify lots of different targets on few samples, it is not necessary to multiplex as you won't save time. Optimizing a multiplex reaction is time-consuming, so it is worth doing it if you plan to work on a large cohort of samples, and always on the same targets.
- If you have strict requirement for your design (cross linked species...), avoid multiplexing, as you can observe annealing between probes or primers, if your design is not optimal.
- However if you manage to have a good PCR efficiency for your singleplex assays, if you don't see any complementarities between your probe and primers, and if you routinely work with the same genes, then it is worth optimizing your multiplex reaction.



HOW DO YOU GET STARTED IN qPCR?

Thermocycler

It seems logical, but you need to get a good thermocycler before you start your qPCR assay. You will most certainly start your qPCR project with the thermocycler available in your lab, or the lab nearby. The majority of the thermocyclers on the market now offer similar characteristics. The table shows the main thermocyclers on the market, and can be considered as a starting point in your decision making process.

The following characteristics should always be taken into consideration:

- **Format:** glass capillaries, plastics tubes, 96-well plates or 384-wells plates. Consider this aspect depending on your requirements, and also based on the reagent availability. A plate format is ideal if you have a high throughput requirement. Glass capillaries can reduce the choice of reagents which you may want to use.
- **Number of detection channels:** if you want to multiplex, choose a thermocycler with a large choice of dyes, in order to have the best combination for your assay. You can find on p. 35 the table showing available dye combinations, based on the type of thermocycler.
- **Software analysis:** Consider the simplicity of the software, the format to export data, statistical analysis performed etc. In addition to qPCR analysis, you can also perform High Resolution Meltcurve analysis with some thermocyclers.
- **Time of run:** the newest thermocyclers enable a very short time of run (30-45 minutes) compared to some of the older ones which take 2 hours.
- **Price and flexibility of the offer.** Some suppliers try to tie up the machine with their corresponding disposable and reagents.

Thermocycler	Supplier	Format	Number of channels	Optical Excitation system	Notes	More information on:
GeneAmp® SDS 5700	Applied Biosystems	96-well plate	1	Halogen lamp	Sales discontinued	http://www.appliedbiosystems.com/
ABI Prism® SDS 7000	Applied Biosystems	96-well plate	4	Tungsten-halogen lamp	Sales discontinued	http://www.appliedbiosystems.com/
ABI Prism® SDS 7700	Applied Biosystems	96-well plate	4	Laser	Sales discontinued	http://www.appliedbiosystems.com/
ABI Prism® SDS 7900 HT	Applied Biosystems	96-well plate or 384-well plate	4	Laser	Robotic plate loading system for high throughput use Normal or FAST block	http://www.appliedbiosystems.com/
ABI Prism® SDS 7300	Applied Biosystems	96-well plate	4	Tungsten-halogen lamp		http://www.appliedbiosystems.com/
ABI Prism® SDS 7500	Applied Biosystems	96-well plate	5	Tungsten-halogen lamp	Normal or FAST block	http://www.appliedbiosystems.com/
StepOne® / StepOnePlus®	Applied Biosystems	48-well / 96-well	3/4	LEDs		http://www.appliedbiosystems.com/
iCycler IQ®	BioRad	96-well plate	4	Halogen lamp	Sales discontinued	http://www.bio-rad.com/
My iQ®	BioRad	96-well plate	1	Tungsten- halogen lamp	Singleplex FAM/SYBR	http://www.bio-rad.com/
iQ5®	BioRad	96-well plate	5	Tungsten-halogen lamp		http://www.bio-rad.com/
CFX 96®	BioRad	96-well plate	5	LEDs	PCR cyclers converted into Real-Time cyclers	http://www.bio-rad.com/
MiniOpticon®	BioRad	48-well plate	2	LEDs		http://www.bio-rad.com/
DNA Engine Opticon® 1	BioRad	96-well plate	1	LEDs	Sales discontinued	http://www.bio-rad.com/
DNA Engine Opticon® 2	BioRad	96-well plate	2	LEDs		http://www.bio-rad.com/
Chromo 4®	BioRad	96-well plate	4	LEDs		http://www.bio-rad.com/
Mx3000P®	Stratagene	96-well plate	4	Quartz tungsten-halogen lamp	Customizable filters	http://www.stratagene.com/
Mx3005P®	Stratagene	96-well plate	5	Quartz tungsten-halogen lamp	Customizable filters	http://www.stratagene.com/
Mx4000®	Stratagene	96-well plate	4	Quartz tungsten-halogen lamp	Sales discontinued	http://www.stratagene.com/
Mastercycler® ep realplex2	Eppendorf	96-well plate	2	LEDs		http://www.eppendorf.com/
Mastercycler® ep realplex4	Eppendorf	96-well plate	4	LEDs		http://www.eppendorf.com/
LightCycler® 1.5	Roche	32 glass capillaries	3	LEDs		https://www.roche-applied-science.com/
LightCycler® 2	Roche	32 glass capillaries	6	LEDs		https://www.roche-applied-science.com/
LightCycler 480®	Roche	96-well or 384-well plate	6	Xenon lamp	Optional HRM module	https://www.roche-applied-science.com/
Smartcycler®	Cepheid Innovation	Up to 96 plastic tubes.	4	LEDs	Independent block of 16 wells. Up to 6 blocks can be combined.	http://www.cephheid.com/
Rotor-Gene™ 2000 / 3000	Corbett Life Science	16 / 32 or 72 plastics tubes	2/4	LEDs	Sales discontinued	http://www.corbettlifescience.com
Rotor-Gene™ 6000	Corbett Life Science		6	LEDs	Includes HRM module	http://www.corbettlifescience.com

cf. the dye-thermocycler compatibility table on p. 35

Table 2. Main characteristics of existing thermocyclers.



qPCR reagents : formats and components?

When choosing a reagent for your qPCR assay, consider the format and the composition, in order to use one which is adapted to your assay, your machine, and your way of working.

Core kits or MasterMixes?

A Core kit, is a kit that contains all components in separate tubes, so you have to mix them yourself. It gives a maximum of flexibility, as you can optimize the concentration of each component of your assay. A qPCR Core kit contains individual tubes of reaction buffer, dNTP/dUTP mix, $MgCl_2$, *Taq* and eventually SYBR® Green.

On the contrary, a MasterMix is a ready-to-use reagent in which all the components are already mixed in an optimized way. Additional stabilizers are added to enable long-term storage. The concentration of components is perfectly adapted to most of the assays. Consequently, the 2 x mix guarantees a high reproducibility, ease of use, and time saving.

Role of Core kits and MasterMixes components

Slight component changes can significantly impact your assay. It is therefore important to consider them individually, understanding the role of each of them.

Taq Polymerase

A HotStart *Taq* polymerase is inactive at low temperatures (room temperature). Heating at 95 °C for several - usually 5 to 10 - minutes activates the enzyme, and the amplification can begin once the primers are annealed. The enzyme is not active until the entire DNA is denatured.

Two major HotStart modifications exist, the antibody-blocked *Taq* and the chemically-blocked *Taq*. The antibody-blocked *Taq* is inactive because it is bound to a thermolabile inhibitor that is denatured during the initial step of PCR. The chemically-blocked *Taq* provides one clear advantage over the antibody-blocked *Taq*, as it is completely inactive at 60 °C, (the hybridization temperature of primers), thus preventing the formation of non-specific amplification and reducing primer dimer formation.

dNTPs / dUTPs

Some kits contain a blend of dNTPs and dUTPs, other ones contain only dNTPs. Using only dNTPs increases the sensitivity, the reason being that the *Taq* incorporates more easily dNTPs than dUTPs. However, using a mix containing dUTPs brings security to the assay, in case of contamination from a previous PCR product. Thanks to the UNG activity in association with incorporated dUTPs, this contamination can be eliminated.

UNG

The Uracil-N-Glycosylase is an enzyme that hydrolyses all single-stranded and double-stranded DNA containing dUTPs. Consequently, if all PCR amplifications are performed in the presence of a dNTPs/dUTPs blend, by carrying a UNG step before every run it is possible to get rid of any previous PCR product.

ROX

Some thermocyclers require MasterMix containing ROX dye for normalization. This is the case for the ABI and Eppendorf machines, and optional on the Stratagene machines. If you work with such machines, it is easier to work with the ROX dye already incorporated in the MasterMix rather than adding it manually. It guarantees a higher level of reproducibility and homogeneity of your assays.

Fluorescein

For iCycler iQ®, My iQ® and iQ5 machines (BioRad thermocyclers), the normalization method for SYBR® Green assay uses Fluorescein to create a “virtual background”. As in the case for the ROX, it is better and easier to use a MasterMix that contains pre-diluted Fluorescein, guaranteeing higher reproducibility and homogeneity of your assays.

MgCl₂

MgCl₂ is necessary for the Reverse Transcriptase and the *Taq* activity. MgCl₂ concentration in MasterMixes is optimized according to the amount of *Taq* and also the buffer composition. However, it may be necessary sometimes to add MgCl₂ and most MasterMixes include an additional tube of MgCl₂.

Inert colored dye

Some buffers also include an inert colored dye, to enable visualization of the buffer when loading in the wells. This colored dye has no effect on the sensitivity of the assay and is a convenient working tool. Note that such mixes, in combination with white plastic plates, provide better levels of fluorescence and a really easy way of working.



Probe and primer quality

We really cannot insist enough, about the importance of this parameter. Your design is actually the most variable parameter, especially in the case of SYBR® Green assays. It is not uncommon to observe undetectable targets, high backgrounds or bad PCR efficiency, due only to the poor design of the primers and probes. The next chapter details more precisely, the rules for probe and primer design. However, consider also the quality and the purification level of your primers and probes, as this can also influence your results.

The classical way of synthesizing oligonucleotides, is the standard Beta-cyanoethyl synthesis cycle that includes: deblocking the first nucleotide, coupling of the second one, capping of any previous nucleotide that has not been coupled, oxidation to stabilize the growing chain, deblocking, then a new cycle starts. At the end of the oligonucleotide synthesis, the crude product is cleaved from the solid support (CPG or polystyrene beads).

Concerning probe coupling, there are different methods by which a dye can be incorporated into a probe. The choice of the coupling method used depends on the structure of the dye and of the length of the probe.

Manual coupling

Manual coupling can be done via amine or thiol groups, using activated dyes such as Texas Red®, TAMRA, JOE, Rhodamine (R6G) or ROX, BODIPY® and other dyes (Alexa®, Marina Blue®, etc.). These labels can be linked internally to any dT-residue, dR-residue or to either terminus. The fluorescent molecule is linked to the oligonucleotide via a spacer (a C-6 spacer generally) to limit the steric hindrances.

Automatic coupling

Automatic coupling uses the appropriate phosphoramidite to which the dye is already coupled like Cy® dyes, FAM, Fluorescein, HEX, TET and Yakima Yellow®. These labels are coupled during the synthesis at the 5' terminus.

The easiest probes to synthesize and to purify are Double-Dye Oligonucleotide probes, as they are not too long. Moreover, the probes where you can use phosphoramidites (automatic coupling method) are usually the easiest ones to synthesize and purify, as there is only one set-up for the synthesis and one purification step. They also give the highest yields. The probes, which require post-labelling procedures (manual coupling method) are more laborious to synthesize and will generally yield lower amounts of product.

Purification level

Regarding the quality, it is recommended to work with purified primers for qPCR assays (Cartridge Reverse Phase purified or HPLC purified primers are ideal). The first aim of any purification method is to remove the by-products, resulting from the synthesis step (mainly salts). The next goal is to enrich the product in full-length oligonucleotides. Consequently, purified primers are more adapted to the sensitive qPCR tool than simply desalted primers. It is particularly true when working with non-specific detection systems like SYBR® Green.

GUIDELINES FOR PRIMERS AND PROBES DESIGN

General rules

Well-designed primers and probes are a prerequisite for successful qRT-PCR. By using well-designed primers and probes, PCR efficiencies of 100 % can be obtained.

If the following primer and probe design guidelines are taken into account you will achieve high PCR efficiencies, specific PCR products, non co-amplification of gDNA and therefore the most sensitive results.

We do recommend in general using a design software (for example Oligo® Primer Analysis Software) to check for all following criteria.

Most thermocycler softwares now offer tools to help you designing primers with the best characteristics. Some of the best softwares are Beacon Designer, Primer Express, and DNA Star... Some other tools are freely available on the web, for example:

- <http://medgen.ugent.be/rtprimerdb/> (human primer and probe database)
- <http://frontend.bioinfo.rpi.edu/applications/mfold/> (for testing secondary structures)
- <http://www.ebi.ac.uk/~lenov/meltinghome.html> (Tm calculators)
- http://frodo.wi.mit.edu/cgi-bin/primer/primer3_www.cgi (Design Primer, also hybridization probes)
- <http://bibiserv.techfak.uni-bielefeld.de/genefisher2> (Design Primer including degenerate primers)
- <http://www.premierbiosoft.com/qpcr/index>

Primers or SNP databases are also freely available, for example OMIM (from the NCBI), SNPedia, Huge Navigator and RT Primer DB (Pattyn F. *et al.*, 2006)



Design guidelines for SYBR® Green I assays

As SYBR® Green I binds to any dsDNA it is important to avoid primer-dimers and non-specific products in SYBR® Green I assays. This can only be avoided by carefully selecting primers that only bind to the selected target. By selecting amplicons between 100 and 150 bp a high level of fluorescence can be obtained without compromising the PCR efficiency.

One important characteristic of primers and probe is the T_m (melting temperature), the temperature at which 50 % of the oligonucleotide is hybridized. The simplified calculation method for the T_m is $T_m = 2(\text{number A+T}) + 4(\text{number G+C})$.

Primers

- Length : 18-30 bases
- GC content : 30-80 % (ideally 40-60 %)
- T_m : 55-60 °C, ΔT_m difference between forward primer and reverse primer should be ≤ 4 °C
- Avoid mismatches between primers and target, especially towards the 3' end of the primer. Avoid runs of identical nucleotides, especially of 3 or more Gs or Cs at the 3' end.
- Avoid 3' end T (allows mismatching)
- Avoid complementarities within the primers to avoid hairpins
- Avoid complementarities between the primers to avoid primer dimers, especially of 2 or more bases at the 3' ends of the primers
- Check if primers are unique and specific (check with a BLAST search: www.ncbi.nlm.nih.gov/BLAST/)
- Design intron spanning or intron flanking primers to prevent amplification of contaminating genomic DNA.

For intron spanning primers the first half of the oligo must hybridize to the 3' end of one exon and to the 5' end of the other exon. In this way only cDNA will be amplified and not gDNA (fig. 11).

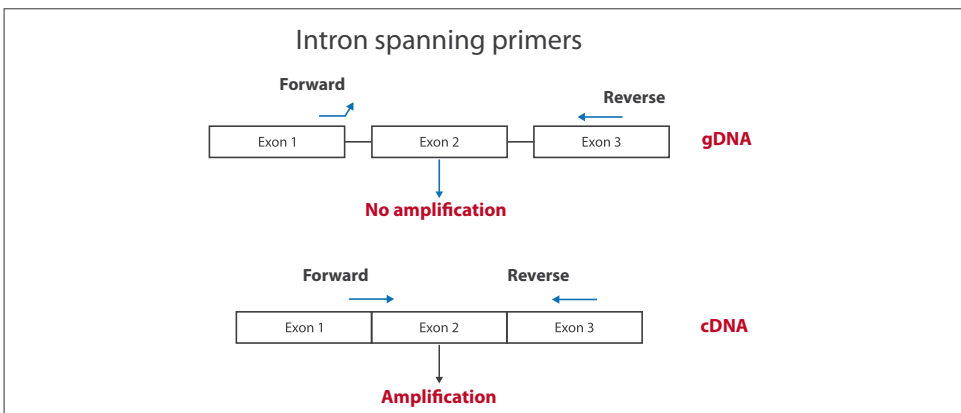


Figure 11. Intron spanning primers

For intron flanking primers, the forward primer must hybridize to one exon and the reverse primer to the other exon. Amplicons from cDNA, without introns, will be smaller than the amplicons from gDNA, which will contain the intron. The bigger amplicon will be amplified less efficiently. The difference in size of these amplicons can be determined via meltcurve analysis. If genes of only one exon are studied contamination with gDNA can only be avoided by DNase treatment of the RNA with RNase free DNase (Vandesompele, 2002).

A recent study indicates that primers with a 5' flap improve Real-Time PCR results (Afonina, 2007). A unique validated 12-bases long AT-rich sequence added to the 5' extremity of both primers might increase the level of fluorescence and improve the sensitivity of the assay. The addition of this sequence is claimed to be especially useful for short primers (when the design is constrained) or on damaged DNA (for example sodium bisulfite treated DNA). The positive effect has been observed both for SYBR® Green I and Double-Dye probe assays.

Amplicon

- Length : 80-150 bp
Shorter amplicons will give higher PCR efficiencies
Longer amplicons will give a higher ΔR_n as more SYBR® green I is incorporated
- GC content
30-80 % (ideally 40-60 %)
- Avoid secondary structures in the amplicon (check with Mfold: <http://frontend.bioinfo.rpi.edu/applications/mfold/>).



Design guidelines for Double-Dye probe assays

In probe assays, primer dimers and non-specific products will not be detected, however, they will influence the PCR dynamics and efficiency.

Therefore, in probe assays they should be avoided as much as possible. For probe assays the amplicons should be kept as short as possible, with the 5' end of the probe as close as possible to the 3' end of the forward primer in case the probe is on the same strand, and as close as possible to the 3' end of the reverse primer in case the probe is on the opposite strand. In this way the 5' nuclease reaction will be optimal.

Experience has showed that it is easier to first design the probe, and then the primers, rather than the other way around.

Primers and probes

Primers

Design rules for a probe assay are the same than for the SYBR® Green assay, see above.

Probes

- Length 18-30 bases, optimal length: 20
Lengths over 30 bases are possible, but it is recommended to position the quencher not at the 3' end, but internally 18-25 bases from the 5' end (normally coupled to a T)
- GC content : 30-80 %
- T_m of the probe must be 8-10 °C higher than the T_m of the primers (8 °C for genotyping, 10 °C for expression profiling)
- Select the strand that gives the probe more Cs than Gs
- Place the probe as close as possible to the primers without overlapping them
- Avoid mismatch between probe and target
- Avoid complementarities within the probe
- Avoid runs of identical nucleotides, especially of 4 or more Gs
- Avoid 5' end G (as this can quench several fluorophores, including FAM)

For multiplex assays / genotyping

- Position the polymorphism in the center of the probe
- Adjust the probe length so that both probes have the same T_m

Amplicon

- Length
80-120 bp optimal
Shorter amplicons will give higher PCR efficiencies and more efficient 5' nuclease reactions
- GC content : 30-80 % (ideally 40-60 %)
Avoid secondary structures in the amplicon (check *with* Mfold:
<http://frontend.bioinfo.rpi.edu/applications/mfold>).

Quenching principle

Fluorophore and quencher: definitions

Many commonly employed techniques for the detection of nucleic acid sequences in a homogenous manner use fluorescence as the signaling technology. Typically a single-stranded probe is labelled with a fluorophore and a quencher molecule. Changes in quenching of the fluorophore, caused by hybridization of the probe to its target nucleic acid, lead to signal generation.

Thus, Real-Time PCR relies not only on the choice of one technology, but also in the choice of the right fluorophore-quencher pairs, whether the assay would be singleplex or multiplex. The choice of a fluorophore and its combination with a quencher, will give different results in terms of sensitivity of the assay.

- A fluorophore is a molecule that emits light of a certain wavelength after having first absorbed light of a specific but shorter wavelength. The emission wavelength is always higher than the absorption wavelength.
- A quencher is a molecule that accepts energy from a fluorophore in the form of light and dissipates this energy either in the form of light or heat.

Different ways of quenching

Quenchers are molecules that can accept energy from a fluorophore and dissipate it by two mechanisms, called proximal, and FRET quenching. A fluorophore absorbs light energy and is promoted to an excited state. In the absence of a quencher the fluorophore falls back to the ground state and the excess of energy is released as fluorescence.

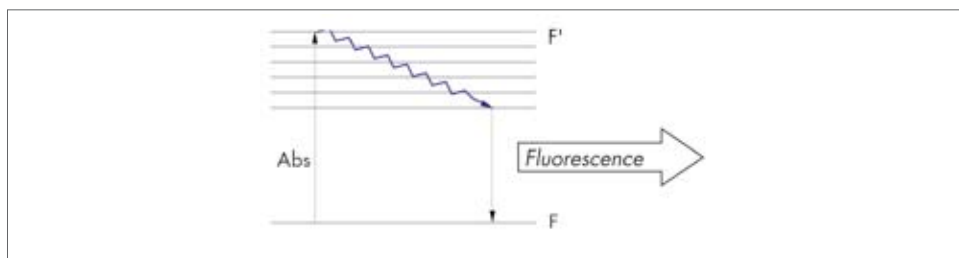


Figure 12. Fluorescence principle

Proximal quenching

When the fluorophore is in close proximity of a quencher molecule, the energy is transferred from the fluorophore to the quencher, which then dissipates the energy as heat (no fluorescence is observed). It is also known as collisional quenching. This type of quenching is used in Molecular Beacons and Scorpions® primers.

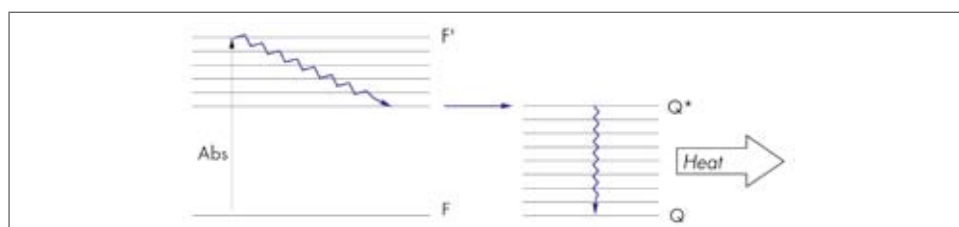


Figure 13. Proximal quenching principle



FRET quenching

The fluorophore transfers its energy to the quencher (which may be another fluorophore); the energy is released from the quencher as fluorescence of a higher wavelength. The efficiency of this process is dependent on the distance between fluorophore and quencher, more precisely on the Förster distance $1/r^6$ (where “r,” is the fluorophore – quencher distance). This type of quenching is used in Double-Dyes probes and LC Hybridization probes. Dark quenchers dissipate the energy as heat, resulting in an absence of fluorescence emitted by the quencher. This type of quenchers allows a better signal-to-noise ratio in the assay.

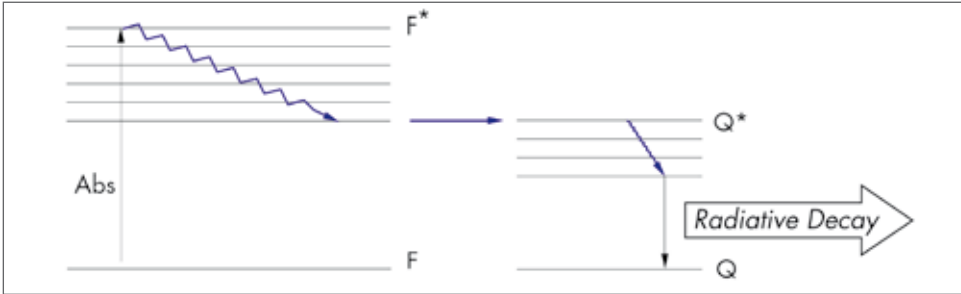


Figure 14. FRET quenching principle

Optimal fluorophore quencher combination

The quencher has to fit the fluorophore

One fluorophore can be combined with multiple different quenchers, but the absorption spectrum of the quencher needs to have good overlap with the emission spectrum of the fluorophore in order to achieve optimal quenching.

- For example, Deep Dark Quencher II absorbs over a large range of the visible spectrum and therefore efficiently quenches most of the commonly used fluorophores, especially those emitting at higher wavelengths.
- Deep Dark Quencher I and Eclipse® Dark Quencher, quench the lower wavelength dyes such as FAM but they are not good at quenching those that emit at high wavelengths such as Cy® dyes.
- The Black Hole Quencher family covers a large range of wavelengths (over the entire visible spectrum and into the near-IR) and is consequently one of the best options for multiplex application.

Quenchers have a quenching capacity throughout their absorption spectrum, but the performance is best close to the absorption maximum.

Quencher	Absorption (nm)	Emission (nm)
TAMRA	565	580
BHQ-1™	534	none
BHQ-2™	580	none
BHQ-3™	670	none
Deep Dark Quencher 1	430	none
Deep Dark Quencher 1	630	none
Dabcyl	475	none

Table 3: Absorbance and emission wavelengths of common Quenchers (BHQ: Black Hole Quencher)

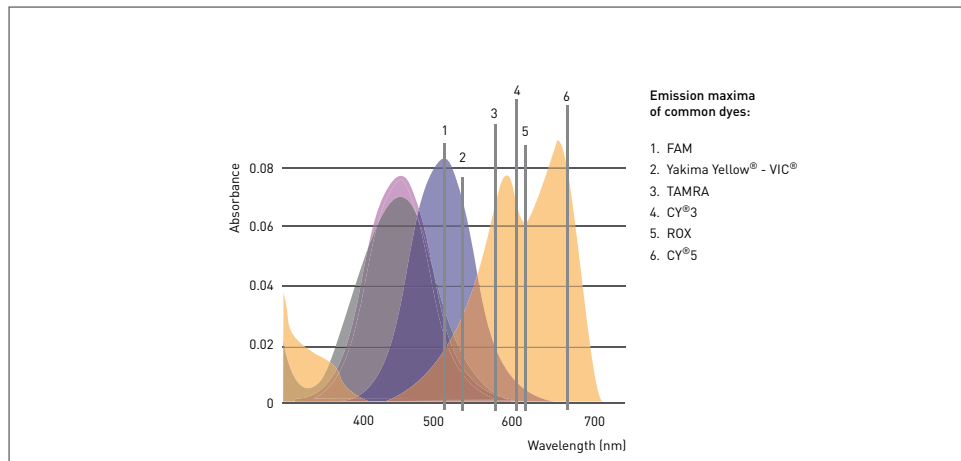


Figure 15. UV-vis absorbance spectra of Deep Dark Quencher (pink), DabCYL (grey), Eclipse® Dark Quencher (purple) and Deep Dark Quencher II (orange)

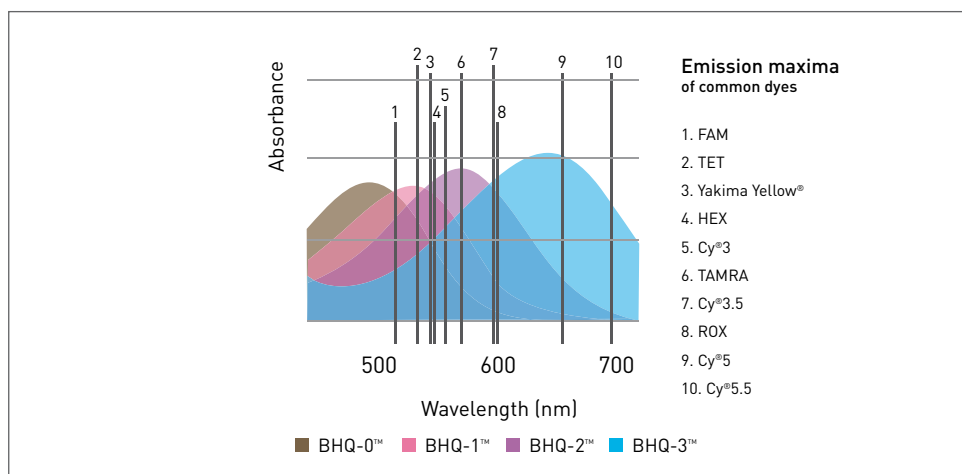


Figure 16. UV absorbance spectra of the Black Hole Quencher™ family

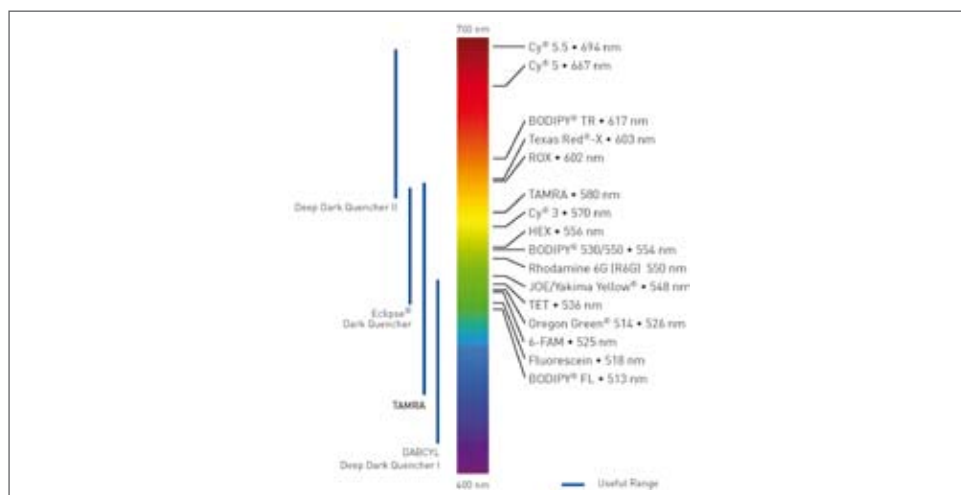


Figure 17. Absorbance range of Dark Quencher dyes and emission maxima for Reporter Dyes



Dye-thermocyclers compatibility table

Almost all detection systems work on almost all Real-Time thermocyclers. The detection channels determine which fluorophores have to be used, and the combination of excitation and detection channels determine which probe systems can be used. It is important to consider, which dyes are the best for multiplexing on the particular machine being used and also any requirements for spectral calibration. This is required because dyes exhibit some spectral overlap, making it necessary to calibrate the machine, which is able to distinguish between the dyes. Spectral calibration kits can be bought for the majority of machines where this may be required.

Double-Dye Oligonucleotides in single-color detection

For single-color detection we recommend to use the combination FAM-TAMRA, FAM-Eclipse® Dark Quencher or FAM-Black Hole Quencher 1. These combinations are the most standard ones, they can be detected on all Real-Time PCR thermocyclers and are easy to synthesize.

Double-Dye Oligonucleotides for multiplexing detection

For multiplex Real-Time PCR, it is very important to select dyes that give a good spectral separation to avoid overlap of the signal. For the first color we recommend the use of FAM-Black Hole Quencher 1. This is the best choice for duplexing or multiplexing PCR, as BHQ-1™ gives a lower background compared to TAMRA.

For the second color, we recommend the use of Yakima Yellow® - BHQ-1™. Yakima Yellow® gives a good spectral separation from FAM, so is the best choice to be combined to FAM-BHQ-1™. Furthermore, it is a good cost effective alternative to VIC, and it can be detected in the same channel at the same wavelength, with no additional calibration needed.

The choice of the third color is dependent on the thermocycler (Table 4). We recommend using either DFO (Dragon Fly Orange), or Texas Red®. These two dyes have both a good spectral separation from FAM and Yakima Yellow®, and both give efficient synthesis with high yields. The choice of the fourth dye is also dependent on the Real-Time thermocycler. We usually recommend the use of Cy® 5, as Cy® 5 has a good spectral separation from FAM, Yakima Yellow® and Texas Red®. Furthermore, this last one also has an efficient synthesis and gives high yields.

Thermocycler	Channel 1	Channel 2	Channel 3	Channel 4	Channel 5	Channel 6	Channel 7
GeneAmp® SDS 5700	FAM						
ABI Prism® SDS 7000	FAM	YYVIC/JOE	DFO/NED	ROX	-		
ABI Prism® SDS 7700	FAM	YYVIC/TET/JOE	DFO/NED	ROX	-		
ABI Prism® SDS 7900 HT	FAM	YYVIC/TET/JOE	DFO/NED	ROX	-		
ABI Prism® SDS 7300	FAM	YYVIC/JOE	DFO/NED	ROX	-		
ABI Prism® SDS 7500	FAM	YYVIC/JOE	DFO/NED/Cy [®] 3	ROX/TR	Cy [®] 5	-	
StepOne®	FAM	YYVIC/JOE	ROX	-	-		
StepOnePlus®	FAM	YYVIC/JOE	DFO/NED	ROX	-		
Cycler iQ®	FAM	YYVIC/HEX/TET/ Cy [®] 3	DFO/NED	ROX	Cy [®] 5		
My iQ®	FAM	-					
iQ5	FAM	YYVIC/HEX/TET	DFO/NED	ROX/TR	Cy [®] 5	-	
CFX 96®	FAM	YYVIC/HEX	ROX / TR	Cy [®] 5	Cy [®] 5.5 / Quasar 705		
MiniOpticon®	FAM	YYVIC/HEX/TET	-	-	-		
DNA Engine Opticon® 1	FAM	-					
DNA Engine Opticon® 2	FAM	YYVIC/HEX/TET	-	-	-		
Chromo 4®	FAM	YYVIC/JOE/TET	DFO/NED/ROX/TR	Cy [®] 5	-		
Mx3000P® (choice of 4 filters)	FAM	TET	YYVIC/JOE/HEX	Cy [®] 3	DFO/NED	TR/ROX	Cy [®] 5
Mx3005P® (choice of 5 filters)	FAM	TET	YYVIC/JOE/HEX	Cy [®] 3	DFO/NED	TR/ROX	Cy [®] 5
Mx4000® (choice of 4 filters)	FAM	TET	YYVIC/JOE/HEX	Cy [®] 3	DFO/NED	TR/ROX	Cy [®] 5
Mastecycler®ep realplex2	FAM	YYVIC/JOE/HEX/TET	-	-	-		
Mastecycler®ep realplex4	FAM	YYVIC/JOE/HEX/TET	DFO/TAMRA	ROX	-		
LightCycler® 1.5	FAM	YYVIC/JOE/HEX/TET	LC Red 705/ Cy [®] 5	-	-		
LightCycler® 2	FAM	LC Red 640/ROX	LC Red 610	LC Red 640	LC Red 670	LC Red 705	
LightCycler 480®	FAM	YYVIC/HEX	LC Red 610	LC Red 640/ROX	Cy [®] 5	LC Cyan 500	
SmartCycler®1	FAM	YYVIC/TET/JOE/ Cy [®] 3	TR	Cy [®] 5			
SmartCycler®2	FAM	Cy [®] 3 / TET	TR/ROX	Cy [®] 5			
Rotor-Gene™ 2000 / 3000	FAM	YYVIC/JOE/TET	ROX/Cy [®] 3.5	Cy [®] 5			
Rotor-Gene™ 6000	FAM	YYVIC/JOE/TET/HEX	ROX/TR/Cy [®] 3.5	Cy [®] 5	Alexa®680/Atto680	Alexa® Fluor-350	

YY: Yakima Yellow® TR: Texas Red® DFO: DragonFly Orange®

Table 4. Dye-thermocycler compatibility table.



EXAMPLE OF A CLASSICAL ASSAY

Some tips and tricks before starting

- Don't forget to include a negative control, and eventually an external or internal positive control (IPC). Choose your normalization method (absolute or relative quantification), *cf.* p. 39.
- Prepare carefully your plate layout, to avoid any doubt while pipetting. No Template Controls should be placed on the plate in such a position that cross contamination is avoided during the set-up; thus, they should be placed if possible away from the highest DNA concentrations, at the top or bottom of the rows, to avoid spilling over the wells, when pipetting.
- Defrost all reagents on ice, and keep SYBR® Green® mix and probes away from light.
- Prepare the reaction mix ideally in a separate room, different from the room where the DNA samples have been prepared, to avoid any contamination.
- Prepare a reaction mix with a volume 5 to 10 % higher than the volume needed. Actually, some pipetting volumes are very small and small droplets can remain outside the tips, leading to loss of reagents. Working in this way can avoid a lack of mix at the end of the plate preparation!
- Don't wait too long after the plate preparation. It is always better to run the qPCR straight after the preparation. You might need to compare two identical plates/samples with different cycling conditions. In this case only, it is better to prepare plates at the same time, then to store one of them in the fridge at 4 °C away from light. You can keep it like this for up to 10 hours.
- When sealing the PCR plates, it is important to ensure that it has been correctly done. If optical films are used, fingerprints and marks should be avoided on the top of the film. Pay special attention to stick the lids to the side of the plate, to avoid evaporation.
- Pipette your sample in triplicate or at least duplicate.
- The plate can be shaken on a plate shaker and spun, before placing it on the machine. This is not an essential step but will ensure that the reaction components are thoroughly mixed and collected at the bottom of the reaction tube. It is useful to also check the wells for bubbles, as bubbles at the bottom of the well can produce unusual plots on the results (when using machines that read from the top).

Classical protocol and thermal profile

By respecting the tips above, prepare a reaction mix containing primers, probes (if applicable), and all the reagents required for the reaction. It is always better to proceed like this because it minimizes differences across the plate and allows for more accurate pipetting.

Commercial MasterMixes often are 2x MasterMixes, meaning half of your reaction volume is MasterMix, and half is composed of primers, probe and template. The primers final concentration should be between 100 nM and 300 nM for a SYBR® Green I assay, and between 100 nM and 900 nM for a probe assay. Probe concentration should be between 100 nM and 300 nM. Nowadays, usual final volume is 25 µl, however it is often possible to work in 20 µl or less to save reagents costs.

Add the DNA sample (template) on the side of each well (water for the NTC), before the reaction mix as it will be transferred down when adding the reaction mix. As the reaction mix is heavier than the DNA, the DNA will be mixed within the reaction mix.

Then, run a first qPCR with the classical thermal cycles on Figure 18.

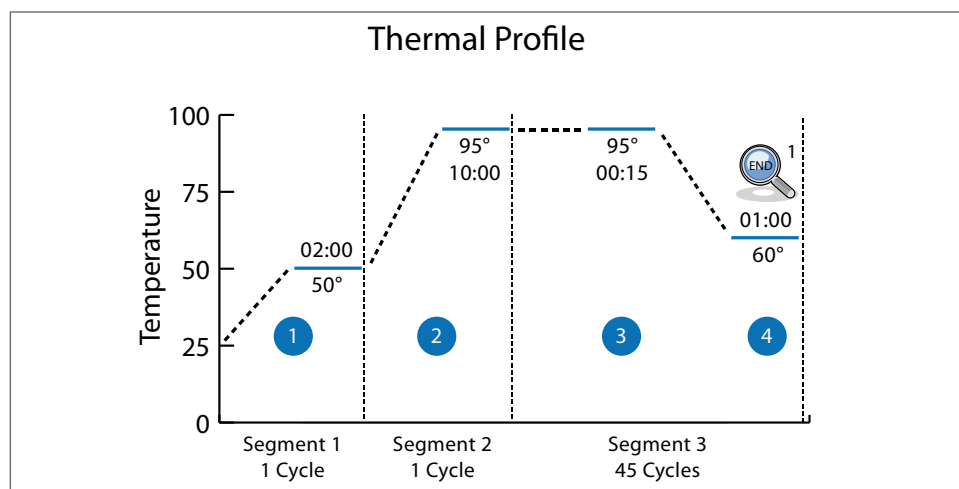


Figure 18. Classical thermal profile including a UNG activation step carry-over prevention step

1. The first step at 50 °C is the UNG step. It is not necessary with all the reagents, as some of them do not include UNG. During this step, the UNG cleaves any contaminating template containing U bases.
2. The 10 minute-step at 95 °C activates the *Taq* and denatures the UNG. Some *Taq* may require shorter activation steps.
3. 15 seconds at 95 °C is the first step of the repeated PCR cycles. The dsDNA template denatures at this temperature.
4. 60 °C for 1 minute, allows the annealing and extension of the primers by the *Taq*.

Depending on your thermocycler, and the characteristics of your primers, the protocol can be shortened and the temperatures adapted. Our optimization tips (p. 49) can help you adapt it to your assay.



Useful definitions

To understand the first results graphs you get, it is necessary to be familiar with some definitions and concepts.

The **Amplification Curve** or **Primary Growth Curve** is usually the first graph we look at. It shows the increase of fluorescence level (called R_n , R_{fu} , $F1\dots$) on the Y axis, compared to the run cycle number on the X axis.

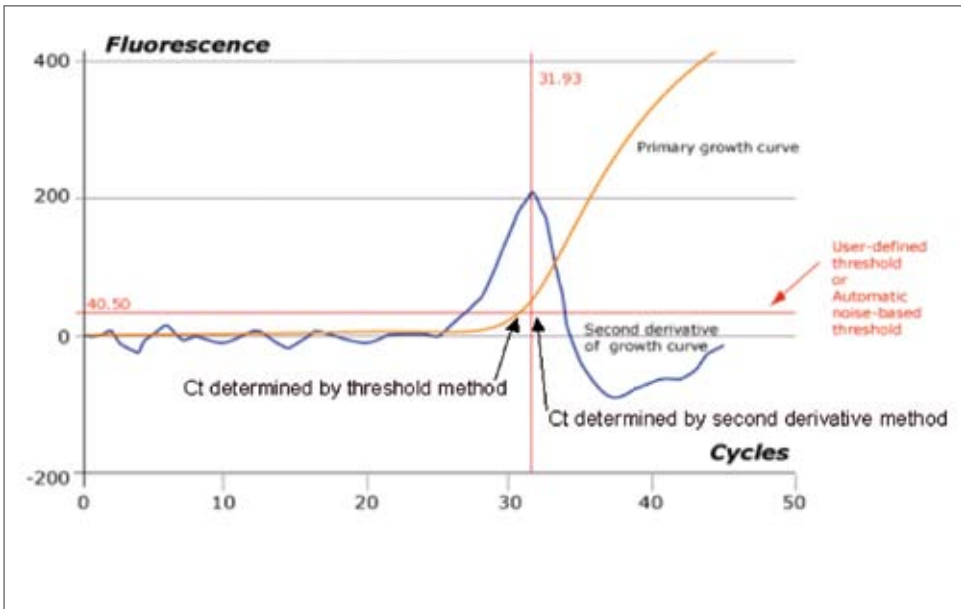


Figure 19. Threshold definition and Ct value calculation

The **baseline** is the average background. It is calculated according to the noise level in the early cycles, when there is no detectable increase in fluorescence, due to PCR products. When the baseline is set up manually, it is important to ensure that only cycles where there is no signal increase are selected !

The **threshold** is the level of fluorescence above the baseline, at which the signal can be considered not to be background. The threshold can be calculated automatically or can be set manually. The automatic calculation of the threshold corresponds to the average baseline + "X" standard deviation of this baseline.

The **Ct value** (C_p for Roche's thermocyclers) is defined as the cycle in which there is a significant increase in reporter signal, above the threshold, i.e. the cycle in which the growth curve crosses the threshold. It is consequently related to the initial amount of DNA and shows also the sensitivity of the assay. The Ct value is consequently in inverse proportion to the expression level of the gene. If the Ct value is low, it means the fluorescence crosses the threshold early, meaning that the amount of target in the sample is high. The **Ct** is therefore dependent on the threshold level, and it is then important to compare the threshold value from one run to another, if no normalization method is used.

The **ROX** is a fluorescent dye that is used as a passive reference for some thermocyclers (Eppendorf, ABI Prism®, Stratagene...). This dye is usually spiked into the MasterMix and the reporter fluorescence is normalized to the ROX signal on ABI and Stratagene Mx machines.

The **Rn** (normalized reporter signal) is calculated by dividing the reporter signal by the ROX signal. This signal is then corrected for pipetting inaccuracies, variation in sample volume, bubbles in sample, plastics inconsistency...

When correction for background variations is required, the ΔRn is used. The ΔRn is calculated by subtracting the normalized baseline from the normalized reporter signal.

$\Delta Rn = Rn - Bn$ where:

$Rn = (\text{reporter signal} / \text{ROX passive reference signal})$

$Bn = (\text{baseline signal} / \text{ROX passive reference})$

Only a minor variation in the ROX passive reference will cause a change in the Rn, Bn, ΔRn , threshold and Ct value. Therefore we recommend using MasterMixes that already contain a ROX passive reference to avoid such variations, which might be interpreted as real differences, but in reality are just experimental artifacts.

Normalization and quantification methods

When analyzing and comparing results of Real-Time qPCR assays many researchers are confronted with several uncontrolled variables, which can lead to misinterpretation of the results. Those uncontrolled variables can be the amount of starting material, enzymatic efficiencies, and differences between tissues, individuals or experimental conditions.

In order to make a good comparison, normalization can be used as a correction method, for these variables.

There are several ways to normalize, which all have some advantages and disadvantages.

The most commonly known and used ways of normalization are : normalization to the original number of cells, normalization to the total RNA mass, normalization to one or more housekeeping genes, and normalization to an internal or external calibrator.

Normalization to number of cells can actually only be done for cell culture and blood samples. In solid tissues and tumors, the amount of cells cannot be counted, only estimated and therefore leads to inaccuracies (Bustin, 2000).

Normalizing to the RNA mass quantity will lead to inaccuracies as the total RNA mass (rRNA, tRNA and mRNA) might contain varying imbalances between rRNA and mRNA (Solinas *et al.*, 2001). The determination of the RNA amount by photospectrometry is another factor causing inaccuracy, as the A260 quantification is strongly influenced by the impurity of the sample, caused by contaminating DNA, free nucleotides and proteins (Bustin, 2000 – Bustin, 2002). As very little is known about total RNA content per cell (in different tissues), a variation between individuals and between normal or highly proliferating tissues, like tumors will cause unwanted inaccuracies (Bustin, 2000 – Bustin, 2002).



Normalization to one housekeeping gene excludes some of the above-mentioned drawbacks, but is also not the perfect method. The advantage is that the variation due to different amounts of RNA can be excluded. However, one must assure oneself, that the housekeeping gene is expressed constantly at the same level throughout the experiment and between samples. It is better then to use several housekeeping genes. A housekeeping gene used for normalization can also be used as a control (*cf.* p. 45).

The two majors methods of normalization are the absolute quantification and the relative quantification (Sellars MJ *et al.*, 2007). Figure 20 summarizes the quantification methods that will be developed further.

In absolute quantification, the exact number of copies of the gene of interest is calculated. In relative quantification, the expression of the gene of interest in a sample is expressed relatively to another gene, another sample, used as a reference.

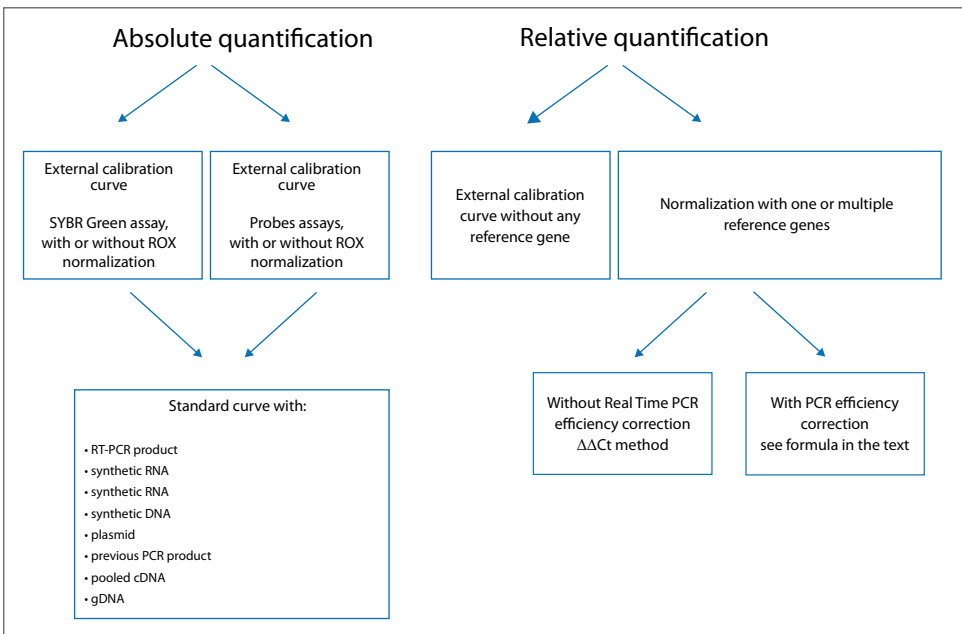


Figure 20: Summary of quantification methods in Real-Time qPCR

Absolute quantification

Absolute quantification requires a standard curve of known copy numbers (*cf.* p. 46 for the construction and analysis of the standard curve). It can be constructed using several standards (Figure 21).

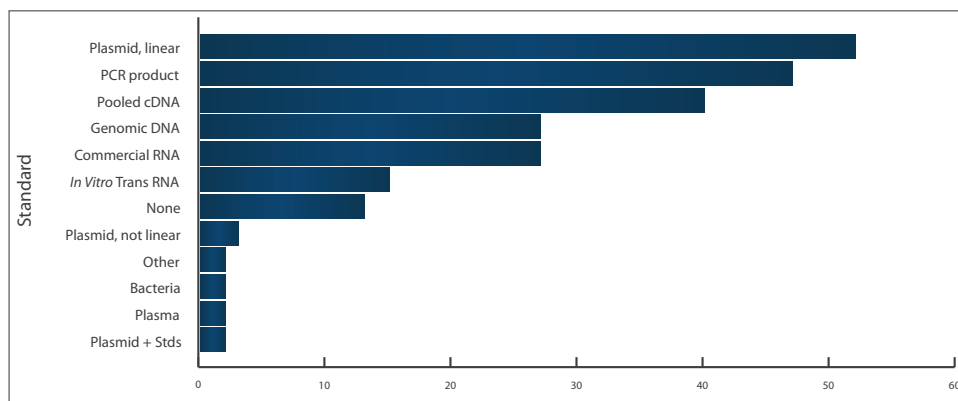


Figure 21: Most frequently used quantification standards. From Nucleic Acid Research Group, (NARG) survey 2007, <http://www.abrf.org/NARG/>

The amplicon being studied can be cloned, or a synthetic oligonucleotide (RNA or DNA) can be used. For gene expression studies, it is recommended to use an RNA standard, so the quantification method takes into consideration any possible RT efficiency variation.

The standard must be amplified using the same primers as the gene of interest and must amplify with the same efficiency. The standards must also be quantified accurately. This can be carried out by reading the absorbance at A260, although this does not distinguish between DNA and RNA, or by using a fluorescent ribonucleic acid stain such as RiboGreen. The problem when using *in vitro* transcribed RNA as standard is that the construction of cDNA plasmid has to be done and has to be reverse transcribed. Moreover, it is not very stable for long-term storage. The advantages of using cDNA plasmid standards is that cDNA plasmids are easy to construct, they can be prepared in large amount and can be stored for a very long time without any difficulties.

Relative quantification

Relative quantification is the most widely used technique. Gene expression levels are calculated by the ratio between the amount of target gene and an endogenous reference gene, which is present in all samples. The reference gene has to be chosen so that its expression does not change under the experimental conditions or between different tissues (Cook NL *et al.*, 2008). There are simple and more complex methods for relative quantification, depending on the PCR efficiency, and the number of reference genes used.



Delta delta Ct ($\Delta\Delta\text{Ct}$) method

This method is the simplest one, as it is a direct comparison of Ct values between the target gene and the reference gene. However, PCR efficiencies of both the target and of the reference gene should be close to 100 % and not differ by more than 10 %. Only the initial experiment requires a standard curve to compare the PCR efficiency of the target and control gene (Schmittgen TD, *et al.*, 2008).

Relative quantification involves the choice of a calibrator sample. The calibrator sample can be the untreated sample, the time=0 sample, or any sample you want to compare your unknown to.

Firstly, the ΔCt between the target gene and the reference gene is calculated for each sample (for the unknown samples and also for the calibrator sample).

$$\Delta\text{Ct} = \text{Ct}_{\text{target}} - \text{Ct}_{\text{reference gene}}$$

Then the difference between the ΔCt of the unknown and the ΔCt of the calibrator is calculated, giving the $\Delta\Delta\text{Ct}$ value:

$$\Delta\Delta\text{Ct} = (\text{Ct}_{\text{target}} - \text{Ct}_{\text{reference}})_{\text{calibrator}} - (\text{Ct}_{\text{target}} - \text{Ct}_{\text{reference}})_{\text{sample}}$$

The normalized target amount in the sample is then equal to $2^{-\Delta\Delta\text{Ct}}$ and this value can be used to compare expression levels in samples.

Different PCR efficiencies

If there is more than 10 % difference in PCR efficiencies, between the reference gene and the target gene, then it is inaccurate to use the $\Delta\Delta\text{Ct}$ method. The value used is then:

Ratio = $(E_{\text{target}})^{\Delta\text{Ct}_{\text{target}}(\text{calibrator} - \text{sample})} / (E_{\text{reference}})^{\Delta\text{Ct}_{\text{reference}}(\text{calibrator} - \text{sample})}$ where :

$$E_{\text{target}} = \text{PCR efficiency of the target}$$

$$\Delta\text{Ct}_{\text{target}(\text{calibrator} - \text{sample})} = (\text{Ct}_{\text{target}})_{\text{in the calibrator}} - (\text{Ct}_{\text{target}})_{\text{in the sample}}$$

When several reference genes are used

As mentioned in the Normalization chapter, it is better to use several control genes, because the expression of one reference gene might change slightly. It is then necessary to use an index reference gene (REF), a method described by M. Pfaffl and J. Vandesompele, 2004.

The normalization by geometrical averaging of multiple internal control genes (Vandesompele *et al.*, 2002) is a robust and innovative approach for accurate normalization. It relies on the use of several housekeeping genes to even out variations in the expression of these genes.

Only the most stable expressed housekeeping genes are taken into the calculation. The calculation can be performed using the freely available software geNorm (<http://medgen.ugent.be/~jvdesomp/genorm/>). The major drawback of this method is the need for many primer pairs and the complicated way to process the data.

EVALUATION OF YOUR qPCR ASSAY

Several parameters can help you validating your assay. As for any experiment, before using a qPCR test routinely, it is important to optimize it, and validate its reproducibility.

Controls

The optimal way of processing is to include several negative and positive controls. The main ones are shown on Figure 22.

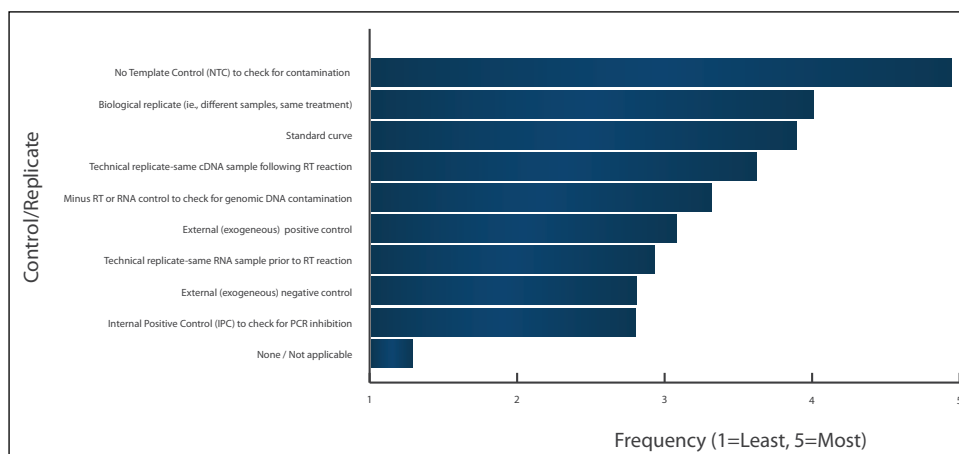


Figure 22: Frequency of replicates controls. From Nucleic Acid Research Group, NARG survey 2007, <http://www.abrf.org/NARG/>

Negative Controls

- **NTC:** the No Template Control contains water and reaction mix, instead of template and reaction mix. Of course you should not observe any amplification in the NTC. However, it is sometimes difficult to achieve. In SYBR® Green assays, if you have many constraints on the design of the primers, you might observe a late amplification due to primer dimers. In any case, the amplification in the non-template control should go up late, and be at least 8 Ct after your specific signal, thus avoiding false positive interpretation. If it is not possible to achieve such a result, it is better to switch to a probe assay that is more specific. A minimum of 2 NTC's should be included.
- **NAC:** the No Amplification Control includes all the components except the Polymerase. It is not possible to do one if you use a MasterMix that contains everything already. However, you can add only probe, primers, template and water, to check the integrity of the probe. If you observe an increase of fluorescence in this case, it is most likely that the probe is degraded. The NAC is less important than the other negative controls, as the probe degradation might also be detected with a high background level.
- **No RT:** the No Reverse Transcriptase control is a sample without the addition of the Reverse transcriptase prior to the PCR step. If the No RT control shows an amplification (and not the NTC control), you can conclude that contaminating genomic DNA is present. The best way to proceed then is to design intron spanning primers or intron flanking primers.



Positive controls

Their main goals are to detect the quality of the reagents used, the presence of inhibitors in your sample, or damaged sample. Positive controls can also be used for normalization.

Endogenous Positive Control

An endogenous control is a second target for which the expression level is measured in the same sample than the Gene of Interest (GOI=target). It can be run in multiplex or in a different well. Endogenous positive control gives indication on the buffer or polymerase action.

An endogenous control gene can be used to standardize the amount of sample RNA or DNA added to the reaction. It should be a stable expressed gene through the experiment, that is used to normalize the results of a variable target gene, and to correct for sample-to-sample variations.

Exogenous Positive Control

It is a different sample in which, the gene of interest is known to be expressed. It can be a previous PCR product, a plasmid containing the gene of interest, *in vitro* RNA transcripts... Thanks to this control, it is possible to check the quality of the primers, probe, and reagents but won't give any indication concerning the sample (degradation, inhibitors...).

Spiking Control or Internal Positive Control (IPC): it is a DNA or RNA that is spiked into the sample and amplified using specific primers. Any variation of this spiking control expression would show that inhibitors are present in the sample (Scipioni A. *et al.*, 2008).

qPCR controls	
➔	<p>NEGATIVE</p> <ul style="list-style-type: none"> • No Template Control • No Amplification Control • No RT Control
	<p>AIM</p> <ul style="list-style-type: none"> Detection of primers dimers and contamination Detection of probe's degradation Detection of genomic DNA contamination
➔	<p>POSITIVE</p> <ul style="list-style-type: none"> • Endogenous Control (same sample, different target) • Exogenous Control (same target, different sample) • Spiking Control (additional DNA spiked into the sample, different target)
	<p>AIM</p> <ul style="list-style-type: none"> Check quality of reagents. Also used for normalization. Check quality of reagents. Detect inhibitors presence Reject false negative in diagnostic assays

Figure 23. Summary of the main existing controls in qPCR and their interest.

Commonly used control kits

A good control should have a constant level of expression between individuals, among different tissues of an organism, at all stages of development, and should not be affected by the experiment treatment. The control gene should also be expressed at similar level as the gene of interest and the range of linear amplification should be known (Bustin, 2000). If a control gene is up and down regulated by the experimental intervention, this will lead to an incorrect normalization and thus to misinterpretation of the results. Therefore, proper validation of presumed stability of expression of the control genes, should be done before studying the target gene.

The table below describes the expression level of commonly used control genes.

Symbol	Accession number	Name	Pseudogene	Relative expression level
-	M10098	18S rRNA0		Very high
-	M11167	28S rRNA0		Very high
ACTB	NM_00101	Beta actin	+	High
GAPD	NM_002046	Glyceraldehyde-3-phosphate dehydrogenase	+	High
UBC	M26880	Ubiquitin C	-	High
B2M	NM_003194	Beta-2-microglobulin	-	High
YWHAZ	NM_000190	Tyrosine 3-monooxygenase activation protein, Zeta polypeptide (Phospholipase A2)	+	Medium
RPL13A	NM_012423	Ribosomal protein L13a	+	Medium
SDHA	NM_004168	Succinate dehydrogenase complex, subunit A	+	Medium
HPRT1	NM_0000194	Hypoxanthine phosphoribosyl-transferase I	+	Medium
TBP	NM_003194	TAT box binding protein	-	Low
HMBS	NM_000190	Hydroxymethyl-bilane synthase		Low

Table 5. Relative expression level and presence (+) or absence (-) of pseudogenes in the genome concerning the most commonly used control genes

Singleplex or multiplex assay for the controls?

The endogenous control can be quantified in a singleplex reaction. The gene of interest and control are, in this case, amplified in separate reactions, using the same cDNA sample. It can also be used in a multiplex reaction, which prevents pipetting errors between the control reaction and gene of interest reaction, which can occur when using singleplex reactions. In this case, normalizing in the same well, will account for the variability in individual PCR assays, like inhibitors. Moreover, multiplex assays have the advantage of being cost effective in terms of reagents required.



Role of the standard curve and evaluation of its parameters

Preparing a standard curve for each gene, which needs to be analysed, can provide a good idea of the performance of the qPCR. The standard curve should cover the complete range of expected expression. Using standard material (plasmid, PCR product etc...), the standard curve should include at least 5 points of dilution, each of them in duplicate (at least). The 10-fold or 2-fold dilution range should cover the largest range of expression levels. Plotting these points on a standard curve, will determine the linearity, the efficiency, the sensitivity and the reproducibility of your assay.

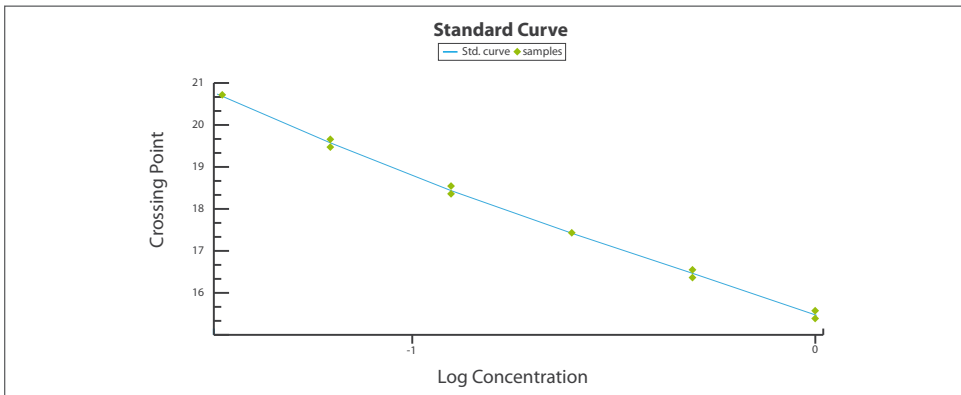


Figure 24: Example of a 2 times-serial diluted cDNA standard curve. Detection using SYBR® Green I on a LC480 thermocycler.

PCR efficiency

The slope of the standard curve gives the efficiency of the PCR reaction by the following equations:

$$\text{Exponential amplification} = 10^{(-1/\text{slope})}$$

$$\text{Efficiency} = 10^{(-1/\text{slope})} - 1$$

If the slope of the standard curve is -3.32 then the PCR is 100 % efficient.

With 100 % efficiency, a 2x dilution gives a ΔCt of 1 between each dilution (each cycle the amount of amplification is doubled).

With a 100 % efficiency, a 10x dilution gives a ΔCt of 3.2 values between each dilution (every 3.2 cycles the amount of amplification is 10 fold higher).

PCR efficiency between 90 %-110 % is acceptable (i.e. a slope between 3.1 and 3.58)

The following factors influence the PCR efficiency:

- Length of the amplicon
- GC content of the amplicon
- Secondary structures in primers, and / or probe, and / or amplicon
- Concentration of reagents

R square (R²)

The R² is a parameter, which indicates how well the data points lie on line. It shows consequently the linearity of the PCR assay. If the R² < 0.95 there is an indication that either the reactions have not been pipetted accurately, or that there is no linear relation between the Ct and the 10 log of the DNA concentration. The latter can be caused by inhibitory factors that are diluted out. A R² higher than 0.985 is acceptable.

Sensitivity

You can compare the sensitivity of your assay by comparing the Ct values. The lower your Ct value is, the higher your assay sensitivity. However, be careful when you compare an assay in terms of sensitivity. You must compare it with the same sample, under the same running conditions, using the same threshold level. A dispersion of the duplicates/triplicates at the lower concentrated sample on the standard curve, can also indicate a lack of sensitivity.

Reproducibility

The reproducibility is indicated by the replicates. The standard deviation between the replicates can be calculated. It's ideally better to avoid replicates that have more than 0.5 Ct of difference

Specificity of a SYBR® Green I assay

For SYBR® Green I assay, an additional point has to be taken into consideration for the evaluation of the assay, this is the specificity. At the end of the run, it's better to systematically run a meltcurve. The temperature is slowly increased from 60 °C to 95 °C whilst continuously monitoring the fluorescence. At a certain temperature, the whole amplified product will fully dissociate, resulting in a drop of fluorescence as the SYBR® Green dissociates from the dsDNA. As the temperature of dissociation is dependant on the length and composition of the amplicon, it is consequently possible to check how many products of amplification are present in the well. A nice meltcurve should show a unique dissociation peak (figure 25). If primer dimers are amplified, a small product is usually observed around 70 °C (figure 26). The primer dimer peak can sometimes even be higher than the peak of the specific product. The primer dimer peak will also be observed even in the NTC control, reinforcing the conclusion that the amplification is non-specific.



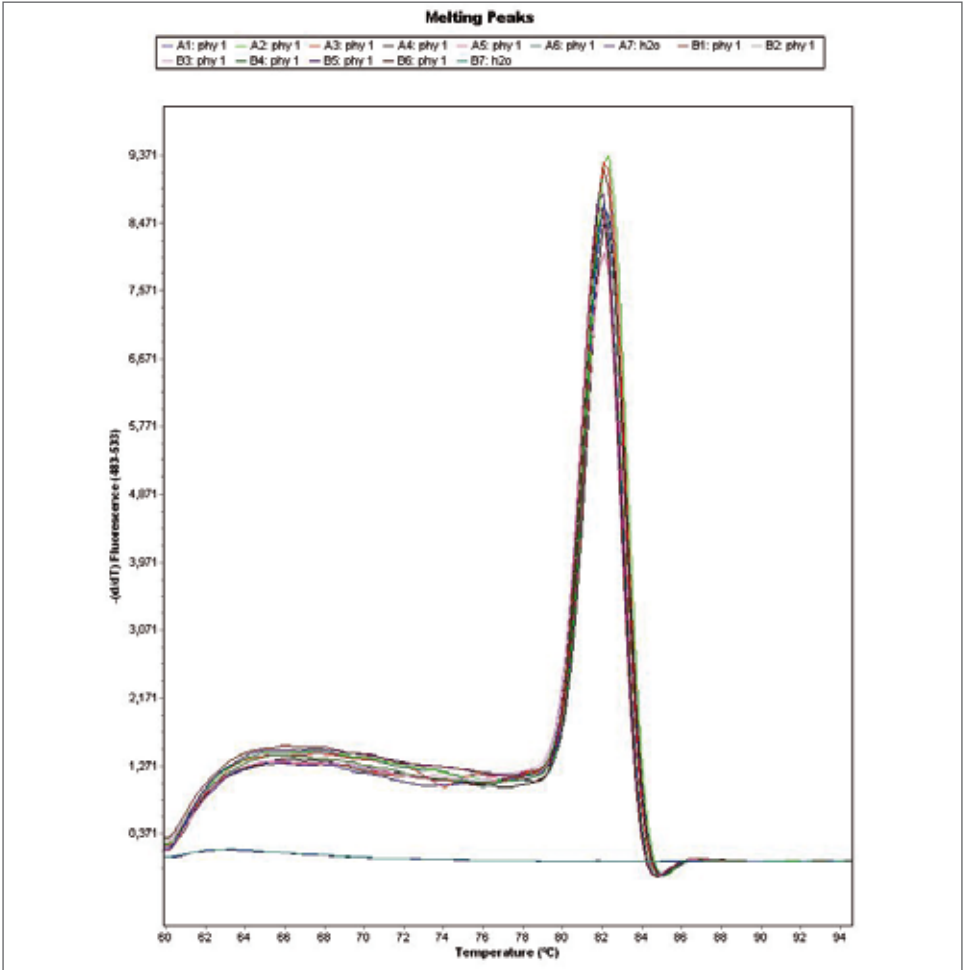


Figure 25. Clean meltcurve (without any primers dimers)

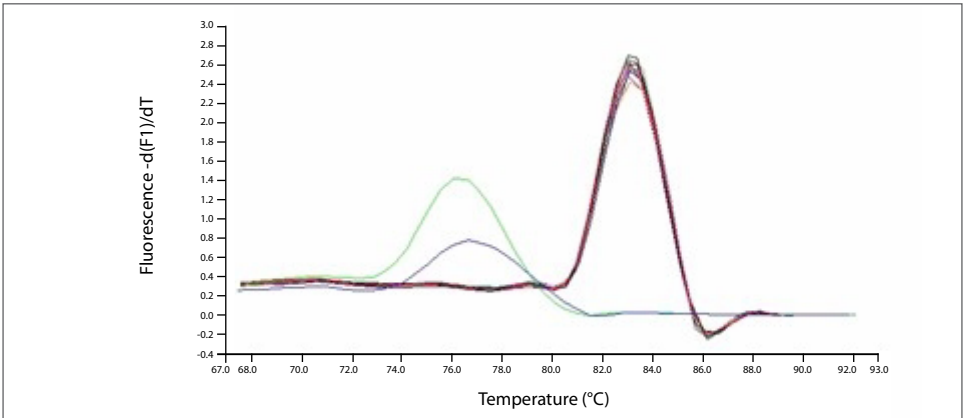


Figure 26. Meltcurve with primers dimers (corresponding to the first peak around 74 °C)

OPTIMIZATION TIPS

In general

If the PCR efficiency is drastically lower than 100 % we would recommend the following:

- Check the primer design and optimize the primer concentration
Trying a new design is sometimes the quickest and cheapest option. To optimize an assay, it is also important to determine the best primer and probe concentration. We recommend testing a matrix of primer and probe concentrations. For SYBR® Green assay: test primers at 50 nM, 100 nM and 300 nM final concentration, each one independently of the other. For probe assay, test primers at 100 nM, 300 nM and 900 nM and test probes at 100 nM and 250 nM. The best combination of forward primer, reverse primer and probe concentrations can then be selected.
- Optimize the temperature reaction and times
Try to change annealing temperature. Start at 60 °C, and then try 58 °C, 56 °C, 54 °C. Try also to increase it up to 64 °C. This annealing temperature of course depends of the T_m of your primers.
- Go for a 3-step instead of a 2-step protocol
In the following cases we recommend performing a 3-step protocol to obtain better results: if results show late Ct values, if the amplicon is quite long, if the results give a less steep growth curve, or if it is not possible to re-do the design.

The protocol will then be as follows:

40 cycles	Denaturation	95 °C for 15 seconds
	Annealing	60 °C for 30 seconds
	Extension	72 °C for 30 seconds

Extension time can be increased with 10-second steps.

It is also possible to add a data collection step after the extension step. To do so, ensure that the detection temperature is at least 3 °C lower than the T_m of the amplicon.



For probe multiplexing assay/qualitative qPCR

General way of processing for a multiplex assay:

- Check carefully the possible complementarity between probes and primers for the different targets.
- Perform (with each target) a singleplex qPCR, with a standard curve to calculate the PCR efficiency.
- Optimize individually each singleplex assays to reach a PCR efficiency as close to 100 %. Multiplex and calculate again the PCR efficiency and the sensitivity of each target. You should not observe a big switch when comparing singleplex and multiplex results for the same target. If there is a loss of sensitivity/efficiency between the singleplex and the multiplex experiment, then try the following things, in this order:
- Limit the probe/primer concentration of the most abundant target.
- Optimize the annealing temperatures by trying steps of 2 °C higher and lower than the current annealing temperature.
- Increase the primer/probe concentration of the less abundant target.

It is actually important to do it in this order, as the main goal is to keep the total primers/probe concentration the lowest possible. A high concentration could partially inhibit the *Taq* polymerase activity.

Such problems can arise when amplification of both genes compete for the PCR reagents. A highly abundant gene can out-compete a less abundant template leading to a bias in the results.

Limiting the primer concentration for the most abundant template can avoid this competition. However, this takes a considerable amount of development. Furthermore, even when working with limited primer concentrations, the detection of rarer target genes may not be possible or may be detectable at much reduced efficiency. Optimizing can also mean adding more units of *Taq* polymerase per reaction and therefore increasing the costs again.

- Be more flexible with your design by using LNA[®] bases. By increasing the T_m of the probe, LNA[®] probes size can be reduced (cf. p. 15) and specificity increased.
- Use a Dark Quencher: when multiplexing, the number of probes in the well increases, and this can lead to an increase of basal fluorescence. Moreover, quenchers such as TAMRA are also fluorophores, leading to high background fluorescence. That is why we recommend the use of Dark Quencher for multiplex assays. There are now several 'Dark quenchers' available. These absorb the energy emitted by the fluorophore, but release it as heat rather than fluorescence, thus leading to a better signal to noise ratio (cf. "Quenching principle"). Dark quenchers are completely non-fluorescent, and have an extremely low background. Using dark quenchers significantly improves signal to noise ratios (higher ΔR_n , a lower background), and thus give a higher sensitivity in multiplex PCR.

TROUBLESHOOTING GUIDE

Contamination

With genomic DNA

The best proof of contamination with genomic DNA is a clean NTC, and a No RT control showing amplification. In this case, carry out a fresh extraction and RT step. If your sample is rare and precious, then use intron spanning primers (*cf.* p. 28).

From a previous PCR run

When using MasterMixes containing a blend of dUTPs and dNTPs, you can easily get rid of a contamination from a previous PCR run. Just use the UNG at the beginning of each run, this will destroy any PCR product containing dUTPs.

Bad PCR efficiency due to primer dimers

Primer dimers in SYBR® Green I assays are easily identifiable on the meltcurve, and the consequence would be a bad PCR efficiency. If you have the choice of the region in which you can design your primers, then the best way to proceed is to do another design. The time wasted for the achieved doubtful result during optimization of an assay with badly designed primers costs much more than a new pair of primers.

If you are limited in the choice of the region for the design, then try to decrease the concentration of primers. You can also try different MgCl₂ concentrations. The best way to proceed however would be to switch to a probe assay which can be more specific.

Inhibition

Decreasing the amount of template can identify the presence of PCR inhibitors in the template. If the Ct values decrease when decreasing amounts of template are used (inverted result), then, it is usually characteristic of the presence of inhibitors. Using an external spiking control can also show the presence of inhibitors.

How to get rid of inhibitors?

There are now several commercial kits to purify RNA samples on columns. You can also try an ethanol precipitation.

If your target is expressed in a sufficient amount, the best way of processing is to dilute the cDNA sample.



Strong background

A strong background can be due to the degradation of the probe. Check it on the NTC, and/or on the multicomponent view, to see if the fluorescence increases without target.

A strong background can also be due to the quencher noise, as quencher such TAMRA is also a fluorophore. Use Dark Quencher, especially for multiplex assays.

A high concentration of DNA can also lead to a strong background, if you do multiplex, try to decrease the concentrations of individual primers and probes.

A strong background can also be explained by simple reasons such as fingerprints on the plate, dirt in the thermocycler block, particles in the tube etc.

Genotyping troubleshooting

Probe cross hybridization or spectral overlap means that, for example, a X probe cross-hybridizes to the Y amplicon. The following solutions can be used:

- Reduce the concentration of the “contaminating” probe, start at 300 nM and decrease to 50 nM.
- Optimize the annealing temperature. A probe that hybridizes despite the SNP, might not anneal at the optimal temperature. Carry out the same test by decreasing and increasing the annealing temperature by 2 °C-steps.
- Increase the concentration of the probe that gives the lowest signal.
- Use LNA[®] to have a more specific X probe, for example by surrounding the SNP by three LNA[®] bases.

It can also be observed that the VIC/Yakima Yellow probe signal overlaps the FAM one, or the contrary. The same solution mentioned above can be tried. You might also try to reverse the fluorophore keeping the same sequences, or to redesign the probes.

FREQUENTLY ASKED QUESTIONS

Is a probe assay more sensitive than a SYBR® Green assay?

Both assays can be equally sensitive, however the probe assay is more specific. Consequently, for difficult to optimize PCR assays, it is recommended to use a probe assay.

How can I get rid of a contamination from a previous PCR product?

Use a kit containing dUTP/dNTP and UNG, which will destroy previous PCR products at the beginning of each PCR program.

Why is my Negative Control, positive?

- Simple contamination by pipetting over the NTC wells: place them as far as possible from the lowest diluted sample (highest concentrated sample).
- The water is contaminated: use freshly purified water
- There are primer dimers in the SYBR® Green assay, leading to a late amplification. This can be checked on the meltcurve.
- The probe is degraded.

How can I avoid gDNA being amplified in a reaction, in which I want to quantify the amount of mRNA?

Perform a DNase treatment of the sample, or design intron spanning primers.

How can I avoid the formation of primer dimers?

Use a HotStart enzyme, optimize the design, and reduce the primer concentration.

Can I use primers I designed for classical PCR for my qPCR experiment?

Some design rules are similar, however one major change is the length of the amplicon, which should ideally not exceed 150 bp in qPCR and can be much longer in PCR.

Why do I have several peaks on my meltcurve?

A peak before the specific peak is usually primer dimers. Two distinct melting peaks can be due to co-amplification of genomic DNA and cDNA, or can be due to a non specific product.

What amount of template should I use in each well?

50 ng of template per well is usually enough. The maximum recommended amount of cDNA is 500 ng in 50 µl, as a too high quantity of RNA/cDNA might inhibit the qPCR reaction.



How can I set up the threshold for my qPCR assay?

The threshold, as the baseline, is set up automatically on most machines, and calculated during the early cycles of the qPCR run. However, you can also set it up manually.

What is an acceptable PCR efficiency for my assay?

A PCR efficiency between 96 % and 104 % is ideal. However a PCR efficiency between 90 % and 110 % is acceptable, and the assay can usually be optimized to achieve this.

What is an acceptable R² value for my assay?

A R² higher than 0.985 is acceptable.

Should I use a kit containing ROX or not?

It depends on your machine. The ROX passive reference is mainly used on the ABI Prism® and Eppendorf machines to normalize the signal. Check our thermocycler-kit compatibility table to be sure you use the right kit for your machine type.

What should I give Eurogentec to get optimization / troubleshooting advice on my results?

Try to sum up the aim of your qPCR assay and the species on which you work. Then you should send us the plate set up, the thermal cycling conditions, the kit and machine you use, the amplification curve and associated Ct values, the standard curve, the meltcurve for SYBR® Green assay, the raw data and multicomponent view if possible on your machine, the sequence and T_m of your primers/probe, the concentration of primers/probe, and the normalization method. Please send all this informations to qPCR@eurogentec.com.

How should I store kits, primers and probes?

Probe MasterMixes are stable for 2 years at -20 °C, SYBR® Green MasterMixes are stable for 1 year at -20 °C. It is also possible to store them at 4 °C for one month. Avoid doing more than 5 cycles of freeze-thawing. RT mixes are stable for one year at -20 °C or -80 °C.

Probes and primers in solution should also be stored at -20 °C in the dark, and are stable for months.

In which buffer should probes be resuspended?

Probes are very sensitive to hydrolysis, so resuspend them in TE 0.01 or just water. Any acidic solution will hydrolyze the probes and give a lower signal to noise ratio.

In which case is a one-step qRT-PCR recommended?

For high-throughput screening and optimized reactions, and to avoid any contamination, a One-Step experiment is the best way to proceed.

EUROGENTEC PRODUCTS

Kits and consumables

Our MasterMixes are developed to work with standard temperature profiles, to be completely adapted to any assay. However, it is recommended to generate a primer optimization matrix and a primer and probe ratio matrix, as this is a crucial step to obtain the lowest signal-to-noise ratio, the earliest Ct value and to save maximum reagents.

If you don't find the reagent you need in the table below, please contact us at qPCR@eurogentec.com.



qPCR Kits for SYBR® Assays • Compatibility with Real-Time thermocycler									
			ABI GeneAmp® SDS 5700 ABI Prism® SDS 7000 ABI Prism® 7300/7700 ABI Prism® 7900HT StepOne® Plus MasterCycler® ep realplex	ABI Prism® 7500	iCycler IQ® My IQ® iQ5	Mx4000® Mx3000P® Mx3005P®	Mx4000® Mx3000P® Mx3005P®	Rotor-Gene™ 2000/3000/6000 MiniOpticon® DNA Engine Opticon® 1 & 2 Chromo 4® CFX96 SmartCycler® 1 & 2 Quantitect® LC480	Sample reference
qPCR MasterMixes with ROX, passive reference	Reference	#RXNs 2.5 µl							
qPCR MasterMix for SYBR® Green I (5 x 1.5 ml)	RT-SN2X-03T	600	✓						RT-SN2X-005+
qPCR MasterMix Plus for SYBR® Green I dTTP, 7.5 ml	RT-SN2X-03+WOUN	600	✓						(1)
qPCR MasterMix Plus for SYBR® Green I w/o UNG, 7.5 ml	RT-SN2X-03WOU+	600	✓						(1)
qPCR MasterMix Plus 7.5 ml for SYBR® Green I	RT-SN2X-03+	600	✓						RT-SN2X-005+
qPCR MasterMix Plus 15 ml for SYBR® Green I	RT-SN2X-06+	1200	✓						RT-SN2X-005+
qPCR MasterMix Plus 50 ml for SYBR® Green I	RT-SN2X-20+	4000	✓						RT-SN2X-005+
One step qRT-PCR MasterMix for SYBR® Green I, 7.5 ml	RT-SNRT-032X	600	✓						RT-SNRT-0052X
Two step qRT-PCR MasterMix for SYBR® Green I	RT-SN2X-03RT	600	✓						(1)
MESA GREEN qPCR MasterMix Plus for SYBR® Assay dTTP, 7.5 ml	RT-SY2X-03+WOUN	600	✓						RT-SY2X-005+WOUN
MESA GREEN qPCR MasterMix Plus for SYBR® Assay, 7.5 ml	RT-SY2X-03+WOU	600	✓						RT-SY2X-005+WOU
MESA GREEN qPCR MasterMix Plus for SYBR® Assay, 15 ml	RT-SY2X-06+WOU	1200	✓						RT-SY2X-005+WOU
MESA GREEN qPCR MasterMix Plus for SYBR® Assay, 50 ml	RT-SY2X-20+WOU	4000	✓						RT-SY2X-005+WOU
MESA BLUE qPCR MasterMix Plus for SYBR® Assay, 7.5 ml	RT-SY2X-03+WOUN	600	✓						RT-SY2X-005+WOUN
MESA FAST qPCR MasterMix Plus for SYBR® Assay, 7.5 ml	RT-SY2X-03+WOUF	600	✓						RT-SY2X-005+WOUF
One step MESA GREEN qRT-PCR MasterMix for SYBR® Assay, 7.5 ml	RT-SYRT-032X	600	✓						RT-SYRT-0052X

(1) Use the RT-SN2X-005+ sample to get a valid preview of the MasterMixes and Core Kits results • dTTP kit will perform slightly better than this sample kit • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.



qPCR Kits for SYBR® Assays • Compatibility with Real-Time thermocycler									
qPCR Core kits with ROX passive reference	Reference	#RX/Is 25 µl	ABI GeneAmp® SDS 5700 ABI Prism® SDS 7000 ABI Prism® 7300/7700 ABI Prism® 7900HT StepOne® Plus MasterCycler® ep realplex	ABI Prism® 7500 My IQ® iQ5	Mx4000® Mx3000P® Mx3005P®	Mx4000® Mx3000P® Mx3005P®	Rotor-Gene™ 2000/3000/6000 MiniOpticon® DNA Engine Opticon® 1 & 2 Chromo 4® CFX96 SmartCycler® 1 & 2 Quantica® LC480	Sample reference	
qPCR Core kit for SYBR SYBR® Green I	RT-SN10-05	1000	✓					(1)	
qPCR Core kit for SYBR SYBR® Green I dTTP	RT-SN10-05NU	1000	✓					(1)	
Two step qRT-PCR Core kit for SYBR® Green I	RT-SMRT-05	1000	✓			with ROX normalisation	without ROX normalisation	(1)	
qPCR MasterMixes with Low ROX passive reference									
qPCR MasterMix Plus for SYBR® Green I Low ROX, 7.5 ml	RT-SN2X-03+WOU LR	600		✓		✓		/	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay Low ROX, 7.5 ml	RT-SY2X-03+WOU LR	600		✓		✓		RT-SY2X-005+WOU LR	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay Low ROX, 15 ml	RT-SY2X-06+WOU LR	1200		✓		✓		RT-SY2X-005+WOU LR	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay Low ROX, 50 ml	RT-SY2X-20+WOU LR	4000		✓		✓		RT-SY2X-005+WOU LR	
MESA BLUE qPCR MasterMix Plus for SYBR® Assay Low ROX, 7.5 ml	RT-SY2X-03+WOU LR B	600		✓		✓		RT-SY2X-005+WOU LR B	
MESA FAST qPCR MasterMix Plus for SYBR® Assay Low ROX, 7.5 ml	RT-SY2X-03+WOU LR F	600		✓		✓		RT-SY2X-005+WOU LR F	
qPCR MasterMixes without ROX passive reference									
qPCR MasterMix Plus for SYBR® Green I No ROX, 7.5 ml	RT-SN2X-03+NR	600				✓	✓	RT-SN2X-005+NR	
Two step qRT-PCR MasterMix for SYBR® Green I No ROX	RT-SN2X-03NRR T	600				✓	✓	(2)	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay No ROX, 7.5 ml	RT-SY2X-03+NRWOU	600				✓	✓	RT-SY2X-005+NRWOU	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay No ROX, 15 ml	RT-SY2X-06+NRWOU	1200				✓	✓	RT-SY2X-005+NRWOU	
MESA GREEN qPCR MasterMix Plus for SYBR® Assay No ROX, 50 ml	RT-SY2X-20+NRWOU	4000				✓	✓	RT-SY2X-005+NRWOU	

(1) Use the RT-SN2X-005+ sample to get a valid preview of the MasterMixes and Core Kits results • dTTP kit will perform slightly better than this sample kit • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.

(2) Use the RT-SN2X-005+NR sample to get a valid preview of the MasterMixes and Core Kits results • dTTP kit will perform slightly better than this sample kit • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.



qPCR Kits for SYBR® Assays • Compatibility with Real-Time thermocycler

	ABI GeneAmp® SDS 5700 ABI Prism® SDS 7000 ABI Prism® 7300/7700 ABI Prism® 7900HT StepOne® Plus MasterCycler® ep realplex	ABI Prism® 7500	iCycler IQ® My IQ® IQ5	Mx4000® Mx3000P® Mx3005P®	Rotor-Gene™ 2000/3000/6000 MiniOpticon® DNA Engine Opticon® 1 & 2 Chromo 4® CFX96 SmartCycler® 1 & 2 Quantica® LC480	Sample reference
qPCR MasterMixes without ROX passive reference	Reference	#RXNs 2.5 µl		with ROX normalisation	without ROX normalisation	
MESA BLUE qPCR MasterMix Plus for SYBR® Assay No ROX, 7.5 ml	RT-SY2X-03+NRW0UB	600		✓	✓	RT-SY2X-005+NRW0UB
MESA FAST qPCR MasterMix Plus for SYBR® Assay No ROX, 7.5 ml	RT-SY2X-03+NRW0UF	600		✓	✓	RT-SY2X-005+NRW0UF
One step MESA GREEN qRT-PCR MasterMix for SYBR Assay No ROX, 7.5 ml	RT-SYRT-032XNR	600	(a)	(a)	✓	RT-SYRT-0052XNR
qPCR Core kits without ROX passive reference						
qPCR Core kit for SYBR Green I No ROX	RT-SM10-05NR	1000	(a)	(a)	✓	(1)
Two step qRT-PCR Core kit for SYBR® Green I No ROX	RT-SMRT-05NR	1000	(a)	(a)	✓	(1)
qPCR MasterMixes with fluorescein						
qPCR MasterMix Plus for SYBR® Green I with fluorescein, 7.5 ml	RT-SM2X-03+NRFL	600				/
MESA Green qPCR MasterMix Plus for SYBR® Assay w/ fluorescein, 7.5 ml	RT-SY2X-03+W0UFL	600		✓		RT-SY2X-03+W0UFL
MESA Green qPCR MasterMix Plus for SYBR® Assay w/ fluorescein, 15 ml	RT-SY2X-06+W0UFL	1200		✓		RT-SY2X-03+W0UFL
MESA Green qPCR MasterMix Plus for SYBR® Assay w/ fluorescein, 50 ml	RT-SY2X-20+W0UFL	4000		✓		RT-SY2X-03+W0UFL
MESA BLUE qPCR MasterMix Plus for SYBR® Assay w/fluorescein, 7.5 ml	RT-SY2X-03+W0UFLB	600				RT-SY2X-005+W0UFLB

(1) Use the RT-SN2X-005+NR sample to get a valid preview of the MasterMixes and Core Kits results. qTTP kit will perform slightly better than this sample kit. For RT reagents, combine this sample to our RT enzyme by simply adding the RT to the mastermix at the recommended concentration.

(a) In combination with ROX passive reference • RT-PARE-03

(b) In combination with fluorescein additive • RT-FLUO-ADD



qPCR Kits for Probe Assays • Compatibility with Real-Time thermocycler									
qPCR MasterMixes with ROX passive reference	Reference	#RXNs 2.5 µl							Sample reference
qPCR MasterMix, 5 x 15 ml	RT-QP2X-03	600	✓						RT-QP2X-005+
qPCR MasterMix Plus, 7.5 ml	RT-QP2X-03-075+	600	✓						RT-QP2X-005+
qPCR MasterMix Plus, 15 ml	RT-QP2X-03+15+	1200	✓						RT-QP2X-005+
qPCR MasterMix Plus, 50 ml	RT-QP2X-03-50+	4000	✓						RT-QP2X-005+
FAST qPCR MasterMix Plus, 7.5 ml	RT-QP2X-03+F	600	✓						RT-QP2X-005+F
FAST BLUE qPCR MasterMix Plus, 7.5 ml	RT-QP2X-03+FB	600	✓						RT-QP2X-005+FB
qPCR MasterMix Plus w/o UNG, 7.5 ml	RT-QP2X-03W0U+	600	✓						RT-QP2X-005+W0U
qPCR MasterMix Plus, 7.5ml QGS	RT-QP2X-03+QGS	600	✓						RT-QP2X-005+QGS
qPCR MasterMix Plus dTTP, 7.5 ml	RT-QP2X-03+W0UN	600	✓						RT-QP2X-005+W0UN
FAST qPCR MasterMix Plus dTTP, 7.5 ml	RT-QP2X-03+W0UNF	600	✓						RT-QP2X-005+W0UNF
FAST BLUE qPCR MasterMix Plus dTTP, 7.5 ml	RT-QP2X-03+W0UNFB	600	✓						RT-QP2X-005+W0UNFB
One step qRT-PCR MasterMix, 7.5 ml	RT-QPRT-032X	600	✓						RT-QPRT-0052X
Two step qRT-PCR MasterMix	RT-QP2X-03RT	600	✓						(1)
qPCR Core kits with ROX passive reference									
qPCR Core kit	RT-QP73-05	1000	✓						(1)
qPCR Core kit dTTP	RT-QP10-05NU	1000	✓						(1)

(1) Use the RT-QP2X-005+ sample to get a valid preview of the MasterMixes and Core Kits results • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.



qPCR Kits for Probe Assays • Compatibility with Real-Time thermocycler									
				ABI Prism® 7500	iCycler iQ® My iQ® iQ5	Mx4000® Mx3000P® Mx3005P®	Mx4000® Mx3000P® Mx3005P®	Rotor-Gene™ 2000/3000/6000 MiniOpticon® DNA Engine Opticon™ 1 & 2 Chromo 4® CFX96 Smartcycle® 1 & 2 Quantitect® LC480	Sample reference
qPCR Core kits with ROX passive reference									
qPCR Core kit only dUTP									(1)
Two step qRT-PCR Core kit									(1)
qPCR MasterMixes with Low ROX passive reference									
qPCR MasterMix Plus Low ROX, 7.5 ml									RT-QP2X-005+WOUJLR
FAST qPCR MasterMix Plus Low ROX, 7.5 ml									RT-QP2X-005+WOUJRF
FAST BLUE qPCR MasterMix Plus Low ROX, 7.5 ml									RT-QP2X-005+WOUJRFB
One step qRT-PCR MasterMix Low ROX, 7.5 ml									(2)
qPCR MasterMixes without ROX passive reference									
qPCR MasterMix Plus No ROX, 7.5 ml									RT-QP2X-005+NR
FAST qPCR MasterMix Plus No ROX, 7.5 ml									RT-QP2X-005+NRFB
FAST BLUE qPCR MasterMix Plus No ROX, 7.5 ml									RT-QP2X-005+NRFB
One step qRT-PCR MasterMix No ROX, 7.5 ml									RT-QPRT-0052XNR
Two step qRT-PCR MasterMix No ROX									(3)
qPCR Core kits without ROX passive reference									
qPCR Core kit No ROX									(3)
Two step qRT-PCR Core kit No ROX									(3)

(1) Use the RT-QP2X-005+ sample to get a valid preview of the MasterMixes and Core Kits results • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.

(2) Use the RT-QPRT-0052XNR sample and add ROX passive reference to the reaction buffer at the recommended concentration (see ROX TDS).

(3) Use the RT-QP2X-005+NR sample to get a valid preview of the MasterMixes and Core Kits results • For RT reagents, combine this sample to our RT enzyme (RT-0125-ER) by simply adding the RT to the MasterMix sample at the recommended concentration.

(a) In combination with ROX passive reference • RT-PARE-03

qPCR miscellaneous					
Fluorescein additive	RT-FLUO-ADD	1 ml	UNG (Uracil-N-Glycosylase)	RT-0610-03	300 U
ROX passive reference	RT-PARE-03	300 µl	UNG (Uracil-N-Glycosylase)	RT-0610-15	1500 U
Reverse Transcriptase Core kit	RT-RTCK-03	600 RXNs (25 µl)	EuroScript/RNase inhibitor mix	RT-0125-ER	3750 U
Reverse Transcriptase Core kit	RT-RTCK-05	1000 RXNs (25 µl)			
qPCR Internal Positive Control					
Internal Positive Control Yakima Yellow ** - TAMRA			RT-PCY-T02		200 RXNs (50 µl)
Internal Positive Control Yakima Yellow ** - BHQ1			RT-PCY-B02		200 RXNs (50 µl)

* For all One-Step RT kits

** Equivalent to VIC®

qPCR validated internal controls												
Target gene	18S rRNA	28S rRNA	β-actin	GAPDH	Ubiquitin C	β-2 microglobulin	Phospholipase A2	HMBS	RPL13a	SDHA	HPRT1	TBP
Reference	RT-CKFT-18S RT-CKYD-18S	RT-CKYD-28S	RT-CKYD-ACTB	RT-CKYD-GAPD	RT-CKYD-UBC	RT-CKYD-B2M	RT-CKYD-YWHAZ	RT-CKYD-HMBS	RT-CKYD-RPL13A	RT-CKYD-SDHA	RT-CKYD-HPRT1	RT-CKYD-TBP
Relative expression level	Very high	Very high	High	High	High	High	Medium	Medium	Medium	Medium	Low	Low
cDNA	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
gDNA	✓	✓	✓	✓	✓	✓			✓			
Organism	H, Ra, M, Rb	H, M	H	H	H	H	H	H	H	H	H	H
Amplicon length	121 bp	87 bp	91 bp	86 bp	76 bp	112 bp	93 bp	119 bp	92 bp	118 bp	100 bp	145 bp

H : human • Ra : rat • M : mouse • Rb : rabbit
 CKFT : FAM-TAMRA probe
 CKYD : Yakima Yellow (VIC® equivalent) - Eclipse® Dark Quencher probe

Each kit contains enough reagents for up to 250 - 50 µl reactions.
 Other dye-quencher combinations are available on request.



96-well plate compatibility with Real-Time thermocyclers		Cover	
LC480			
Quantica®			
MiniOpticon®			
Chromo 4			
DNA Engine Opticon® 2			
DNA Engine Opticon® 1			
MasterCycler ep realplex			
Mx3005P®			
Mx3000P®			
Mx4000®			
iQ5			
My iQ®			
iCycler iQ®			
ABI Step One Plus			›
ABI Step One			›
ABI Prism® SDS 7900 FAST			›
ABI Prism® SDS 7900HT		›	›
ABI Prism® SDS 7700		›	›
ABI Prism® SDS 7500 FAST			›
ABI Prism® SDS 7500		›	›
ABI Prism® SDS 7300		›	›
ABI Prism® SDS 7000		›	›
ABI GeneAmp® 5700		›	›
qPCR 96-well plate type			
Sub-skirted, high profile, frosted or white type (ABI systems)		Flat caps (RT-FLAT-300)	
RT-PL96-AB		Adhesive seals (RT-OPSL-XX**)	
RT-PL96-ABW			
Sub-skirted, low profile, frosted or white type (ABI FAST systems)		Adhesive seals (RT-OPSL-XX**)	
RT-PL96-AF			
RT-PL96-AFW			
Non-skirted, high profile, natural or white type		Flat caps (RT-FLAT-300)	
RT-PL96-MQ		Adhesive seals (RT-OPSL-XX**)	
RT-PL96-MQW			
Non-skirted, low profile, frosted or white type		Flat caps (RT-FLAT-300)	
RT-PL96-OP		Adhesive seals (RT-OPSL-XX**)	
RT-PL96-OPW			
Non-skirted, low profile, white		Flat caps (RT-FLAT-300)	
RT-PL96-LC480		Adhesive seals (RT-OPSL-XX**)	

qPCR optical semi-dome 8-cap strips are also available (RT-OCAP-300; 300 strips). They can be used on several machines.

* Plate must be cut to size (48 wells) to fit Step One or MiniOpticon

** XX: RT-OPSL-25 for 25 seals / RT-OPSL-100 for 100 seals



Primers and probes

Since 1987, Eurogentec has been committed to providing oligonucleotides of the highest quality matching all research needs. Our oligonucleotide synthesis service offers custom synthesis in large variety of scales with a wide range of modifications (FAM, DragonFly Orange®, Yakima Yellow®, TAMRA, BHQ, Deep Dark Quencher...) and purification levels:

- Unmodified and Modified Primer Pairs
- Double-Dye Probes
- LNA® Double-Dye Probes
- Molecular Beacons
- Standard and Duplex Scorpions®
- LC Hybridization Probes
- Plexor Primer Pairs



Eurogentec know-how currently meets the satisfaction of more than 35 000 users worldwide.

And if you became the next Eurogentec's customer for PCR and qPCR oligonucleotides...

If you need help with the design of your oligonucleotides, information about prices and delivery estimation or any other question, please contact our Customer HelpDesk, staffed by qualified and experienced scientists in Molecular Biology and Chemistry at info@eurogentec.com.



REFERENCES

- ABOU-BACAR A., PFAFF A.W., GEORGES S., LETSHER-BRU V., FILISETTI D., VILLARD O., ANTONI E., KLEIN J.P., CONDOLFI E., "Role of NK cells and gamma interferon in transplacental passage of toxoplasma gondii in mouse model of primary infection", *Infection and Immunity*, Vol.72 (3), 1397-1401, 2004
- AFONINA I., ANKOUDINOVA I., MILLS A., LOKHOV S., HUYNH P., MAHONEY W., "Primers with 5' flaps improve real-time PCR", *Biotechniques*, Vol.43, 770-774, 2007.
- BUSTIN S., "Absolute quantification of mRNA using Real-Time reverse transcriptin polymerase chain reaction assays", *Journal of Molecular Endocrinology*, Vol.25, p.169-193, 2000
- BUSTIN S., "Quantification of mRNA using Real-Time reverse transcription PCR (RT-PCR); trends and problems", *Journal of Molecular Endocrinology*, 29, 23-39, 2002
- COOK NL., VINK R., DONKIN JJ., VAN DEN HEUVEL C., "Validation of reference genes for normalization of real-time quantitative RT-PCR data in traumatic brain injury.", *J Neurosci Res.*, 2008.
- ERALI M., VOELKERDING KV., WITTEWIT CT., "High resolution melting applications for clinical laboratory medicine.", *Exp Mol Pathol.* ,85(1):50-8, 2008.
- FLEIGE S. and PFAFFL MW., "RNA integrity and the effect on the real-time qRT-PCR performance", *Mol Aspects Med.* 27(2-3): 126-139, 2006
- FREEMAN W.M., WALKER S.J., VRANA K.E., "Quantitative Rt-PCR:pitfalls and potential", *Biotechniques*, Vol.112 (22), 124-125, 1999
- HALFORD W.P., "The essential prerequisites for quantitative RT-PCR: pitfalls and potential", *Nature Biotechnology*, Vol.17 (19), 835, 1999
- HERNANDEZ M., PLA M., ESTEVE T., PRAT S., PUIDOMENECH P. and FERRANDO A., "A specific real-time quantitative PCR detection system for event MON810 in Maize YieldGard® based on the 3'-transgene integration sequence", *Transgenic Research*, Vol.12, 179-189, 2003
- HOLLAND P., ABRAMSON R.D., WATSON R. and GELLAND D.H., "Detection of specific polymerase chain reaction product by utilizing the 5'-3' exonuclease activity of thermos aquatus.", *Proc. Natl. Acad. Sci. USA*, Vol.88, 7276-7280, 1991
- HOOGWIJS D., HOUTHOOFD K., MATTHIJSSSENS F., VANDESOMPELE J., VANFLETEREN JR., "Selection and validation of a set of reliable reference genes for quantitative sod gene expression analysis in *C. elegans*.", *BMC Mol Biol.* 9:9, 2008
- ILVANDERBROUCKE I., VANDESOMPELE J., DEPAEPE A. and MESSIAEN L., "Quantification of splice variants using Real-Time PCR", *Nucleic Acids Research*, Vol.29 (13), 2001
- ISHII T., SOOTOME H., YAMASHITA K., "Practical evaluation of universal conditions for four-plex quantitative PCR.", *Anal Bioanal Chem.* 388(1):271-8, 2007
- IVANOVA A., ROSCH N., "The structure of LNA:DNA hybrids from molecular dynamics simulations: the effect of locked nucleotides.", *J Phys Chem A.*111(38):9307-19, 2007
- JU J., RUAN C., FULLER C.W., GLAZER A.A. and MATHIES R.A., "Fluorescence energy transfer dye-labeled primers for DNA sequencing and Analysis", *Proc. Natl. Acad. Sci. USA.*, Vol.92, 4347-4351, 1995
- KARKARE S., BHATNAGAR D., "Promising nucleic acid analogs and mimics: characteristic features and applications of PNA, LNA, and morpholino.", *Appl Microbiol Biotechnol.* 71(5):575-86, 2006
- KUBISTA M., ANDRADE JM., BENGTSOON M., FOROOTAN A., JONAK J., LIND K., SINDELKA R., SJOBACK R., SJOGREEN B., STROMBOM L., STAHLBERG A., ZORIC N., "The real-time polymerase chain reaction", *Mol Aspects Med.* 27(2-3):95-125, 2006.
- LE CANN P., RANARIJAONA S., MONOHO S., LEGYADER S., FERRE V., "Quantification of human astroviruses in sewage using real-time RT-PCR", *Res. Microbiol.*, Vol.155 (1), 11-15, 2003
- LETERTRE C., PERELLE S., DILASSE F., ARAR K., FACH P., "Evaluation of the performance of LNA and MGB probes in 5' nuclease assays", *Molecular and cellular probes*, Vol.17 (6), 307-11, 2003
- LI Y., ZHOU X., YE D., "Molecular beacons: an optimal multifunctional biological probe.", *Biochem Biophys Res Commun.* 5:373(4):457-61, 2008
- MASON M., PROVERO P., VAIRA A.M. and PAOLO P., ACCOTTO, "Estimating the number of integrations in transformed plants by quantitative real-time PCR", *BMC Biotechnology*, Vol.2 (20), 2002
- MYAKISHEV M., KHRIPIN Y., MU S., DEAN M., "High-Throughput SNP Genotyping by Allele-Specific PCR with Universal Energy-Transfer-Labeled primers", *Genome Research*, Vol.1, 163-169, 2001
- MULLIS K. and FOBONA F., *In methods in Enzymology*, Vol.155, 335, 1987
- OVERBERGH L., GIULIETTI A., VALCKS D., DECALLONNE B., BOUILLON R. and MATHIEU C. J., "The use of real-time transcriptase PCR for the quantification of Cytokine Gene Expression", *Biomol.Tech.* Vol.14, 33-43, 2003
- PATTYN F., ROBBRECHT P., DE PAEPE F., SPELEMAN F., VANDESOMPELE J., « RTPrimersDB : the real-time PCR primer and probe database, major update 2006". *Nucleic Acids Res* 34, D684-8, 2006

SCHMITTGEN TD., LEE E.J., JIANG J., SARKAR A., YANG L., ELTON T.S., CHEN C., "Real-time PCR quantification of precursor and mature microRNA." *Methods* 44, 31-8, 2008

SCHMITTGEN TD., LIVAK KJ., "Analyzing real-time PCR data by the comparative C(T) method." *Nat Protoc.* 3(6):1101-8, 2008

SCIPIONI A., MAUROY A., ZIANT D., SAEGERMAN C., THIRY E., "A SYBR Green RT-PCR assay in single tube to detect human and bovine noroviruses and control for inhibition.", *Virology* 5(1):94, 2008

SELLARS MJ., VUOCOLO T., LEETON LA., COMAN GJ., DEGNAN BM., PRESTON NP., "Real-time RT-PCR quantification of Kuruma shrimp transcripts: a comparison of relative and absolute quantification procedures.", *J Biotechnol.* 129(3):391-9, 2007

SOLINAS A., BROWN L.J., MXKEEN C., MELLOR J.M., NICOL J., THELWELL N., BROWN T., "Duplex Scorpions primers in SNP analysis and FRET applications", *Nucleic Acids Res.* Vol.29, E96, 2001

TAKACS T., JENEY C., KOVACS L., MOZES J., BENCZIK M., SEBE A., "Molecular beacon-based real-time PCR method for detection of 15 high-risk and 5 low-risk HPV types.", *J Virol Methods.*149(1):153-162, 2008

TYAGI S., MARRAS S.A.E. and KRAMER F.R., "Wavelength-Shifting Molecular Beacons", *Nature Biotechnology*, Vol.18, 1191-1196, 2000

VANDESOMPELE J., DE PRETER K., PATTYN F., POPPE B., VAN ROT N., DE PAEPE A. and SPELEMAN F., "Accurate normalization of real time quantitative RT-PCR data by geometric averaging of multiple internal control genes", *Genome Biology*, Vol.3 (7), 2002



DISCLAIMERS

The polymerase Chain Reaction (PCR) process is covered by patents owned Hoffmann-La Roche, Inc. use of the PCR process requires a license.

A license under U.S. Patents 4,683,202, 4,683,195 and 4,965,188 or their foreign counterparts, owned by Roche Molecular Systems, Inc. and F. Hoffmann-La Roche Ltd (Roche), has an up-front fee component and a running-royalty component. The purchase price of this product includes limited, non-transferrable rights under the running-royalty component to use only this amount of the product to practice the Polymerase Chain Reaction ("PCR") and related processes described in said patents solely for the research and development activities of the purchaser when this product is used in conjunction with a thermal cycler whose use is covered by the up-front fee component. Rights to the up-front fee component must be obtained by the end user in order to have a complete license. These rights under the up-front fee component may be purchased from Applied Biosystems or obtained by purchasing an Authorized Thermal Cycler. No right to perform or offer commercial services of any kind using PCR, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is hereby granted by implication or estoppel. Further information on purchasing licenses to practice the PCR Process may be obtained by contacting the Director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404 or the Licensing Department at Roche Molecular Systems, Inc., 1145 Atlantic Avenue, Alameda, California 94501.

Use of UDG employs U.S. Patents 5,035,996, 5,945,313, 5,683,896 and their foreign counterparts licensed to Eurogentec S.A. from Invitrogen Corporation.

Cyclic-substituted unsymmetrical cyanine dyes are covered by U.S. Patents 5,436,134 and 5,658,751 and licensed to Eurogentec S.A. by Molecular Probes, Inc. in the direct research field.

The use of certain fluorogenic probes in 5' nuclease assays may be covered by U.S. Patents 5,210,015 and 5,487,972, owned by Roche Molecular Systems, Inc., and by U.S. Patent 5,538,848, owned by The Perkin-Elmer Corporation. The purchase of Eurogentec's product does not provide a license to use this patented technology. A license must be obtained by contacting the Director of Licensing Applied Biosystems, 850 Lincoln Centre Drive, Foster City, CA 94404 or the Licensing Department at Roche Molecular Systems Inc., 1145 Atlantic Avenue, Alameda, CA 94501.

The use of Scorpions® Probes for research purposes is covered by a license to Eurogentec S.A. from DxS Ltd. No rights to perform or offer diagnostic, commercial testing or other commercial services for money or money's worth are granted by the supply of this product. Should you wish to use this product for any other purpose not covered by this license, please contact DxS Limited at Manchester Incubator Building, 48 Grafton Street, Manchester, M13 9XX.

The use of Custom Molecular Beacon Probes employs the following patent rights licensed to Eurogentec S.A. from The Public Health Research Institute of the City of New York: [a] the claims of U.S. Patent Application Serial No. 08/439,819; [b] the claims of U.S. Patent Application Serial No. 08/990,176; and [c] all patent application claims and all patent claims, U.S. and foreign, that cover improvements in the foregoing claimed probes, kits and assays and whose subjects are inventions and discoveries that are made by The Public Health Research Institute of New York. The use of Custom Molecular Beacon Probes is limited to the field of use comprising the internal use by an end user of a product solely in assays of the end user [or in applications of the end user's customer, if the end user is performing contract research] in scientific research and development.

The oligonucleotides probes and primers using the Yakima Yellow®, amidite and the Eclipse®, Non-Fluorescent Quencher are for research purposes only, and may not be used for commercial, clinical, diagnostic or any other use. These products or portions thereof are subject to proprietary rights of Epoch Biosciences, Inc. and are made and sold under license from Epoch Biosciences, Inc. under the patent rights WO 0142505 and other pending patent applications. There is no implied license for commercial use with respect to these products. A license must be obtained directly from Epoch Biosciences, Inc. with respect to any proposed commercial use of these products.

LNA oligonucleotides produced under license from Exiqon® A/S for research use only. Not for resale or for therapeutic use or use in humans.

Alexa Fluor® dyes are provided under license from Molecular Probes, Inc., for research use only, and are covered by U.S. Patents 5,696,157 and 6,130,101 and their foreign counterparts.

Texas Red®, Rhodamine Red™, Marina Blue®, dyes are provided under license from Molecular Probes, Inc., for research use only, and are covered by U.S. Patents 4,774,339; 5,132,432; 5,227,487; 5,274,113; 5,696,157; 5,798,276; 5,830,912; 5,846,737; 5,955,612; 6,162,931 and their foreign counterparts.

TRADEMARKS

ABI Prism® is a registered trademark of The Perkin-Elmer Corp. Amplifluor® is a registered trademark of Intergen Alexa®, Alexa Fluor®, Marina Blue®, SYBR®, SYBR gold® and Texas Red® are registered trademarks of Molecular Probes, Inc. Cy® and Cy5 Direct are registered trademarks of Amersham Biosciences Corp. DNA Engine Opticon® is a registered trademark of MJ Research, Inc. Eclipse®, MGB®, MGB Eclipse® and Yakima Yellow® are registered trademarks of Epoch Biosciences, Inc. Exiqon® is a registered trademark of Exiqon A/S GeneAmp® and TaqMan® are registered trademarks of Roche Molecular Systems, Inc. GoldStar® is a registered trademark of Eurogentec S.A. HyBeacons® is a registered trademark of LGC. I-Core® and Smartcycler® are registered trademarks of Cepheid, Inc. iCycler iQ® is a registered trademark of Bio-Rad Laboratories, Inc. Mx3000p® and Mx4000® are registered trademarks of Stratagene Primer Express® is a registered trademark of Applied Biosystems Corporation. Oligo® Primer Analysis Software is a registered trademark of Molecular Biology Insights Inc. Quantica® is a registered trademark of Techne, Ltd and Techne, Inc. Scorpions® is a registered trademark of DxS Ltd Uniprimer® is a registered trademark of Serologicals Royalty Company, Rotor-Gene™ is a registered trademark of Corbett Life Science.



EUROPEAN OFFICES

BELGIUM

Eurogentec Headquarters

Tel.: + 32 4 372 74 00

Fax: + 32 4 372 75 00

info@eurogentec.com

Eurogentec Biologics Department

Tel.: + 32 4 372 74 00

Fax: + 32 4 365 16 04

UNITED KINGDOM

Eurogentec Ltd.

Tel.: + 44 1794 511 411

Fax: + 44 1794 522 417

info.uk@eurogentec.com

FRANCE

Eurogentec France s.a.s.u.

Tel.: + 33 2 41 73 33 73

Fax: + 33 2 41 73 10 26

GERMANY

Eurogentec Deutschland GmbH

Tel.: + 49 221 258 94 55

Fax: + 49 221 258 94 54

THE NETHERLANDS

Eurogentec Nederland b.v.

Tel.: + 31 43 352 06 98

Fax: + 31 43 354 19 65

SWITZERLAND

Eurogentec s.a., succursale de Genève

Français :

Tel.: + 32 4 372 74 00

Fax: + 32 4 372 75 00

Deutsch:

Tel.: + 49 221 258 94 55

Fax: + 49 221 258 94 54

US OFFICE

NORTH AMERICA

Eurogentec North America, Inc.

Production Site

Tel.: + 1 858 793 2661

Toll-free: + 1 877 387 6436

Fax: + 1 858 793 2666

info.usa@eurogentec.com

ASIAN OFFICES

JAPAN

Nippon EGT Co., Ltd.

Tel.: + 81 76 411 02 77

Fax: + 81 76 452 03 99

info@n-egt.com

SINGAPORE

AIT Biotech, Pte. Ltd.

Tel.: + 65 6778 9388

+ 65 6778 6822

Fax: + 65 6778 2962



EUROGENTEC s.a. - LIEGE Science Park - 4102 Seraing - BELGIUM
RPM Liège - VAT (BE) 0427.348.346 - Tel.: + 32 (0)4 372 74 00
Fax: + 32 (0)4 372 75 00 - info@eurogentec.com
www.eurogentec.com